



Lympha Press® Product Catalog



Table of Contents

About Lympha Press®.....03

Lympha Press® Pneumatic Compression Therapy System.....04

Lympha Press® World Famous Garments.....06

Lympha Press® Lower Body Systems.....08

Lympha Press® Upper Body Systems.....10

12 Chamber Compression Therapy Systems.....12

Lympha Press® Optimal Plus.....14

Lympha Press® Mini Press.....15

Lympha Pants®.....16

LymphaPod®.....18

FL Adjustable Leg Garment20

ComfySleeve™ 12 Chamber Leg Garment21

Lympha Jacket™.....22

ComfySleeve™ 1-7524

12-Chamber Arm Garment26

4/8 Chamber Compression Therapy Systems28

Lympha Press® PCD-51™.....30

PCD Pants Lite™.....32

8-Chamber Leg Garment34

4-Chamber Leg Garment35

4-Chamber Arm Garment36

ComfyLite™.....38

Garment Fitting Tips.....40

Measuring Forms.....42

Indications and Contraindications.....46

Quick Reference Guide48

About Lympha Press®

Lympha Press® is the global leader in treatment for lymphedema and venous insufficiency



Expertise

Have been serving 100,000's of users worldwide for over 40 years



Innovation

The world's most innovative line of garments and appliances, with exclusive treatment modes



Research

Over 30 clinical trials and peer-reviewed publications



Recognition

Ranked #1 by patients as the most effective home care treatment for lymphedema



Lympha Press® Dynamic Compression Therapy Systems



Lympha Press® Optimal Plus (Model 912)



Lympha Press® Mini Press (Model 960)



PCD-51™ (Model 731A)

Unique Features

Lympha Press® has the most versatile and advanced dynamic compression systems on the market, featuring patented built-in capabilities and unique treatment modes.

Proven success

Our compression systems have been improving patients' lives around the globe for over 40 year. They are recommended by leading practitioners and therapists.

Advanced Technology


Lympha Press® combines patented technology with groundbreaking research that visually confirms the efficacy of pneumatic compression therapy.

"We have actually seen in real time the way the Lympha Press device can activate the lymphatics to enhance lymphatic flow using near-infrared imaging technology."

"This validates what we understand clinically; that patients improve when they use the Lympha Press® pneumatic compression device (PCD)."

- Caroline E. FIFE, MD, Chief Medical Officer Intellicure, Inc.


Lympha Press® World-Famous Garments

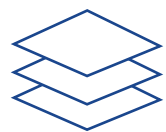


Patented overlapping
chamber design



Up to 24 chambers
for ultimate precision

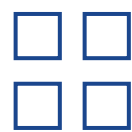

Built-in features
for optimal comfort

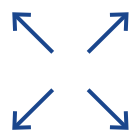

Extended
treatment area



Lympha Press® Lower Body Systems

Leg garments

- Treat leg(s) only
- Used when no torso treatment is desired or needed
- Fits most sizes (expanders are available)

Lympha Pants®

- Treat legs, abdomen, trunk, genitals, hips, and buttocks
- Adjustable boot design fits most leg lengths; expanders are available.

LymphaPod®

- Fits larger patients
- Treats legs, abdomen, trunk, genitals, hips, and buttocks
- Recommended for sizes as shown below:

Waist: Above 200 cm

Hips: Above 180 cm

Upper thigh: Above 95 cm

System compatibility for lower body garments

System	4-Chamber Leg Garment	8-Chamber Leg Garment	ComfySleeve™	FL 12-Chamber Leg Garment	PCD Pants Lite™	Lympha Pants®	LymphaPod®
Optimal Plus (Model 912)			✓	✓		✓	✓
Mini Press (Model 960)			✓	✓		✓	
PCD-51™ (Model 731A)	✓	✓			✓		

Lympha Press® Upper Body Systems

Arm garments

- Treat arm and shoulder
- Used when no torso treatment is desired or needed
- Fits most sizes

Comfysleeve™ 1-75 and Comfylite™


- Treats one arm, axilla, shoulder, back, chest, and abdomen
- Applies compression to one side of torso
- Applies compression to armpit area
- Fits most sizes

Lympha Jacket™

- Treats both arms and torso including shoulder, back, chest, and abdomen
- Applies compression to both sides of torso

System compatibility for upper body garments

System	4-Chamber Arm Garment	12-Chamber Arm Garment	ComfySleeve™ 1-75	ComfyLite™	Lymph Jacket™
Optimal Plus (Model 912)		✓	✓		✓
Mini Press (Model 960)		✓	✓		✓
PCD-51™ (Model 731A)	✓			✓	

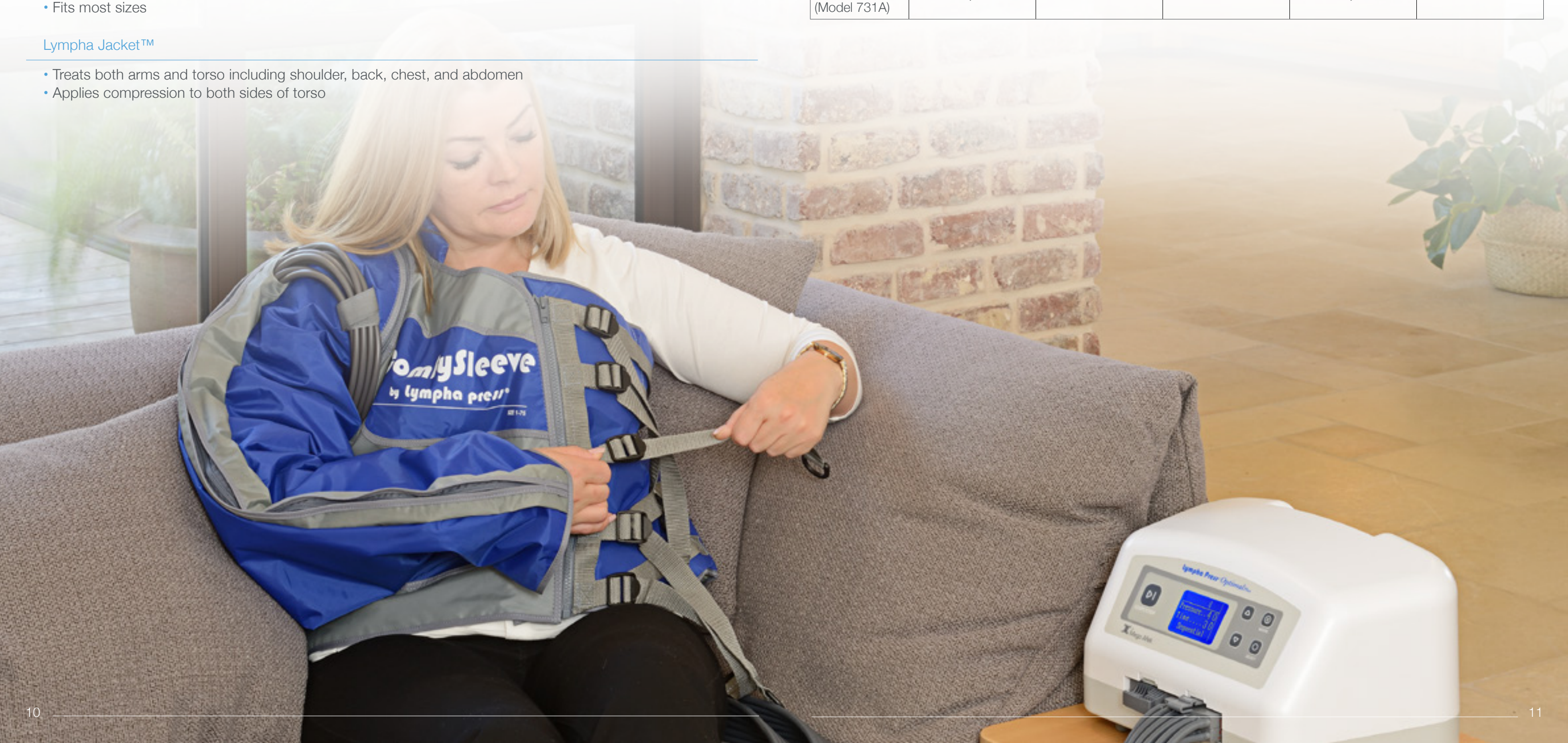

12-Chamber Compression Therapy Systems



Lympha Press® Optimal Plus



Overview

Lympha Press® Optimal Plus, model 912 is the most versatile and advanced dynamic compression system for lymphedema treatment on the market today. It offers unique, state of the art features and treatment modes including:

- Sequential mode for effective reduction of edema.
- Wave mode for comfortable, effective treatment in patients with distal discomfort.
- Pretherapy™ prepares the proximal treated area to receive lymph fluid before the main treatment.
- Post-therapy™ additional therapy to a specific area after the main treatment is completed.

Lympha Press® Optimal Plus includes Bluetooth module for communication with the Lympha Press® Mobile Application. The Lympha Press® Mobile Application can be used with Android and iOS Smartphones for communication with the Optimal Plus device.

Features:

- Adjustable range of 20 to 90 mmHg in a gradient, to enhance distal to proximal flow.
- Automatic adjustable pause time to ensure vein refill
- User-friendly, clear interface
- Programmable via computer, for pressure adjustment of individual chambers
- Compatible with comprehensive range of unique and innovative garments
- Air deflation after treatment
- Full compliance meter. Easy data download via USB port.
- Multi Voltage
- Cost Effective

Technical specifications:

Cycle	Sequential, Wave, Pretherapy, Post-therapy
Pressure Range	20-90 mmHg
Size	34 cm × 24 cm × 19.8 cm (13.3" × 9.4" × 7.8")
Weight	4.5 kg (9.92 lbs.)
Electrical	100-240 VAC, 50-60 Hz

Ordering information

Catalog	Product Name
L10008C2	Lympha Press Optimal Plus 912 with Bluetooth connectivity

For indications and contraindications, see page 48

Lympha Press® Mini Press



Overview

Lympha Press® Mini Press, model 960 offers fully pressure-controlled 12-chamber sequential compression, with user-selected treatment time and pressure (20-80 mmHg). It features basic interface and sturdy, robust design for ease of use and ultimate compliance.

Features:

- Sturdy, reliable and user-friendly for clinic and home use
- Gradient sequential compression
- Multi-voltage
- Compatible with Lympha Press® range of unique and innovative 12-chamber garments

Technical specifications:

Cycle	Sequential
Pressure Range	20-80 mmHg
Size	34 cm × 24 cm × 19 cm
Weight	2.9 kg
Electrical	100-240 VAC, 50-60 Hz

Ordering information

Catalog	Product Name
L1000824	Lympha Press Mini Press 960

For indications and contraindications, see page 48

Overview

The world’s famous Lympha Pants® has extended features for optimal treatment of lymphedema in the lower trunk. Its overlapping 24-chambers compression structure ensures a fluent, comfortable massage and an even pressure application, without gaps.



Features

- Patented overlapping 24-chambers compression structure
- Complete lower body treatment, including abdomen, hips, buttocks and genital area
- For independent home use: the boot structure, once formed, will hold in place
- Patent-protected one-size garment to accommodate most sizes; XL size is available
- Unique Hook-and-Loop fasteners
- Hose bundles exiting either the heel area or the pants’ sides
- Patented “custom boot” feature forms a boot to fit each patient’s leg length
- Expanders are available to increase circumference.
- Zipper opening for catheter access
- Easy to use – as simple as wearing a pair of pants

Ordering Information

Model	Cat. No.	Expanders	Cat. No.
Lympha Pants™ Hip Hose Exit	L4MP2136C	Pants Expanders (pair)	L4Z5LP00
Lympha Pants™ Heel Hose Exit	L4MP2130C		
Lympha Pants™ XL Hip Hose Exit	L4MP2236C		
Lympha Pants™ XL Heel Hose Exit	L4MP2230C		


Compatible with



Lympha Press® Optimal Plus



Lympha Press® Mini Press

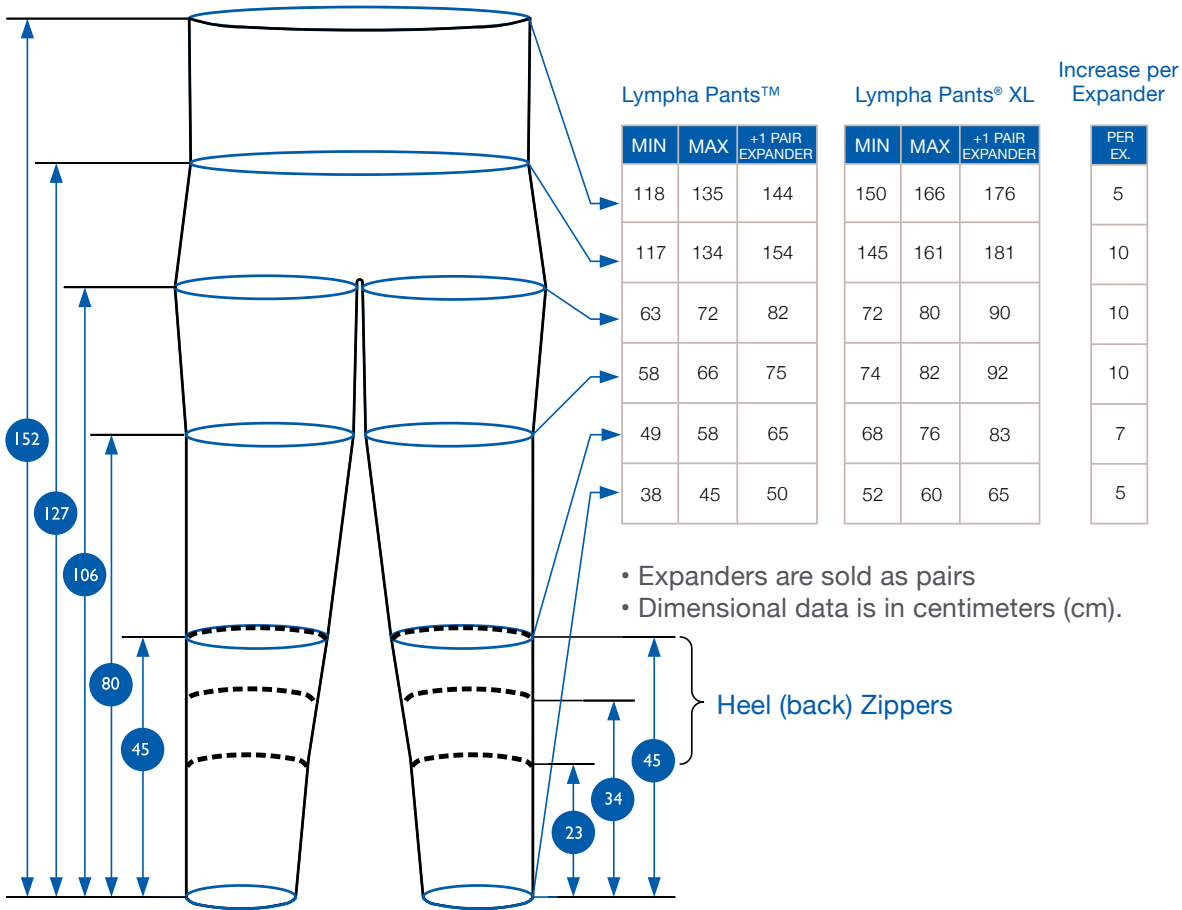

Technical notes

- The garment has three zipper fasteners for circumferential adjustment, with the MIN circumference at the inner zipper position and the MAX circumference at the outer zipper position.



Overview

LymphaPod® is a home care system designed for larger patients. It treats the entire lower body, including legs, hips, and abdomen. It conforms to irregular limbs and shapes and is designed for ease of use by patients with limited mobility.

Features

- Treats the entire lower body including hips, abdomen, and buttocks.
- Special emphasis on the inner thigh to treat lobules.
- Designed for use by patients with limited mobility (some patients may need assistance).
- One size fits most patients from 150 to 330 + Kgs.
- 40 overlapping air chambers.
- Available with hose attachments to the right side or the left side for easy patient access.


Ordering Information

Model	Cat. No.
LymphaPod® Left Hose Exit	L4MP800L
LymphaPod® Right Hose Exit	L4MP800R
LymphaPod® Left Expander	L41580AL
LymphaPod® Right Expander	L41580AR

Compatible with



Lympha Press® Optimal Plus

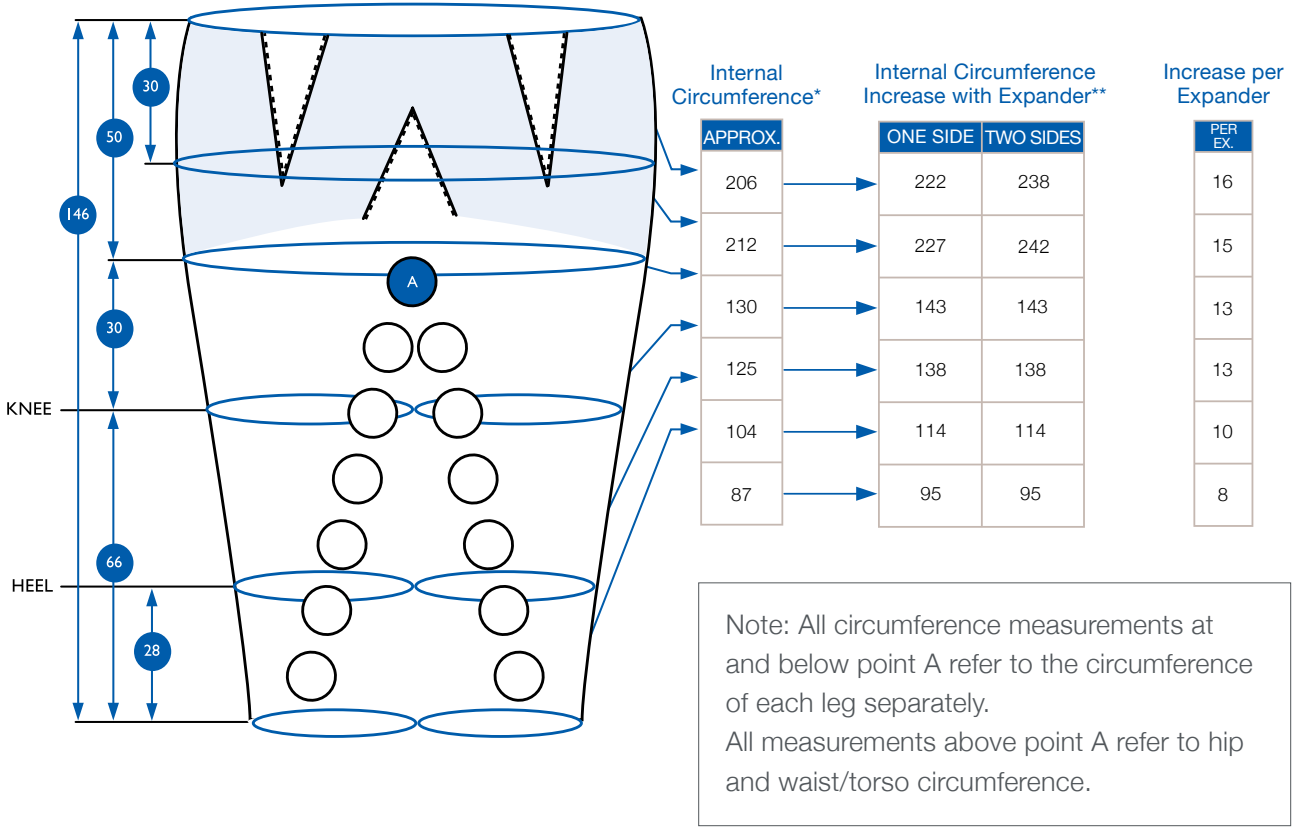

Technical notes

- Lympha Pod® can be used only with Lympha Press® Optimal Plus.
 - Lympha Pod® garments are available with hoses exiting from the right or left side of the garment.
 - Dimensional data is in centimeters (cm).
 - Expanders are sold individually - specify LEFT or RIGHT.
- * Location of patient landmarks will vary depending on height and leg length.
** Leg measurements for side with expander added.



FL Adjustable Leg Garment

Overview

A unique, fully adjustable and conforming garment for treatment of the foot and leg.

Features

- 12 overlapping chambers with four chambers over the foot and ankle, for treatment of distal swelling.
- Fully adjustable with hook-and-loop fasteners for a perfect fit.
- Available in a range of lengths and widths, including extra wide.

Ordering Information

Part No	Model/Size	Width	A - B (Top)	C - D (Ankle)	Length	Instep (E)
L4132075	FL 2-75	Narrow	60 - 70	40 - 50	75	30
L4132085	FL 2-85		66 - 76	40 - 50	85	
L4132095	FL 2-95		70 - 80	40 - 50	95	
L41330105	FL 2-105		74 - 84	40 - 50	105	
L4133075	FL 3-75	Medium	66 - 76	43 - 53	75	34
L4133085	FL 3-85		70 - 80	43 - 53	85	
L4133095	FL 3-95		74 - 85	43 - 53	95	
L41330105	FL 3-105		80 - 90	43 - 53	105	
L4134075	FL 4-75	Wide	78 - 90	56 - 68	75	36
L4134085	FL 4-85		88 - 96	56 - 68	85	
L4134095	FL 4-95		90 - 108	56 - 68	95	
L4Z5A030	FL Expander 30X30 cm. Adds up to 20 cm. circumference.					

Compatible with



Lympha Press® Optimal Plus



Lympha Press® Mini Press

ComfySleeve 12 Chamber Leg Garment



12-chamber garments

Overview

Ideal for patients who need a garment that is easy to put on and remove.

Features

- Overlapping chambers – our patented system for smooth and effective therapy without backflow or gaps.
- Designed for ease of use.
- Wide selection of standard sizes.
- Specially designed chambers create a conforming fit with air.
- Amputee and custom garments are available.

Ordering Information

			Dimensions							With Expander					
Part No	Model/Size	Width	A	B	C	D	E	F	Length	Model	A	B	C	D	E
L41L2275	COMFY 2-75	Narrow	64	54	43	40	43	28	75	BT	75	64	51	45	48
L41L2285	COMFY 2-85		68	56	42	38	41	25	85	CT	80	67	51	45	45
L41L2295	COMFY 2-95		77	64	49	40	42	26	95	ET	89	77	59	49	48
L41L2365	COMFY 3-65	Medium	79	68	56	51	51	36	65	BT	89	77	65	57	55
L41L2375	COMFY 3-75		84	71	57	52	52	39	75	CT	93	81	65	59	57
L41L2385	COMFY 3-85		88	74	60	50	51	37	85	DT	98	84	68	57	56
L41L2475	COMFY 4-75	Wide	94	79	64	64	58	40	75	FT	109	93	76	70	63

Compatible with



Lympha Press® Optimal Plus




Lympha Press® Mini Press

Overview

The Lympha Jacket™ treats both arms together with the entire upper body (chest, shoulders, back, and abdomen) according to the principles of manual lymph drainage.

Features

- Complete bilateral torso therapy.
- Treats both arms simultaneously.
- 24 overlapping chambers prevent reflux and promote directional flow to healthy lymphatics.
- Easy to put on and take off, just like a regular jacket.
- Easy to clean for good hygiene and skin care.

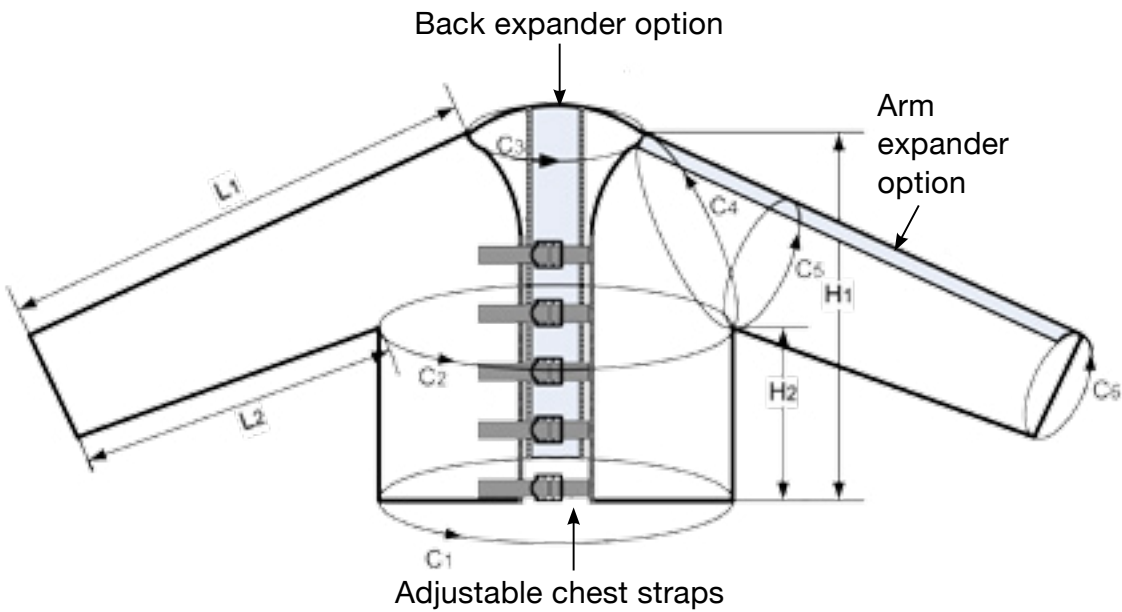

Ordering Information

						Adjustable Front Straps		Arm/Shoulder		
Part No	Model/Size	L1	L2	H1	H2	C1/C2	C3	C4	C5	C6
L4MJ2125	Lympha Jacket, Short Arm	80	57	60	29	105 - 123	48 - 64	68	48	40
L4MJ2135	Lympha Jacket, Long Arm	96	72							

Versatile Expander Options		
L4Z52121	Expander Set for Lympha Jacket, Short Arm. Contains three expanders (1 back expander, 2 arm expanders)	Each expander adds 10 cm to circumference.
L4Z52130	Expander Set for Lympha Jacket, Long Arm. Contains three expanders (1 back expander, 2 arm expanders)	
L4150542	Jacket Expander, Back	
L4150822	Jacket Expander, Short Arm	
L4150823	Jacket Expander, Long Arm	Extra wide back expander for large abdomens
L415T510	Jacket Expander, Triangular	
L4150510	Jacket Expander, Chest	
		Zips on the front, adds 15 cm to torso circumference.

Compatible with



Lympha Press® Optimal Plus



Lympha Press® Mini Press



Overview

ComfySleeve™ 1-75 treats the entire arm, shoulder, and one side of the back, chest, and abdomen.

It provides gentle, targeted therapy to the axilla (armpit), to assist lymphatic drainage.



Features

- Treats arm, chest, shoulder, back, and abdomen.
- Conforms to the patient’s body shape, with targeted compression to the armpit area.
- Easy to put on and remove.
- Easy to clean for good hygiene.
- Modular system is cost effective and easy to use.
- Change the size and side of the garment by switching fasteners.

Ordering Information

MODEL ComfySleeve™ (12 chambers) Select Side Treated	CAT. No. For Left Arm and Torso: L41L0175L Left Arm + Vest For Right Arm and Torso: L41L0175R Right Arm+Vest
Select: Small Fastener	For Left Arm: L42L0FSL For Right Arm: L42L0FSR
Select: Medium Fastener	For Left Arm: L42L0FML For Right Arm: L42L0FMR
Select: Large Fastener	For Left Arm: L42L0FLL For Right Arm: L42L0FLR


Compatible with



Lympha Press® Optimal Plus



Lympha Press® Mini Press

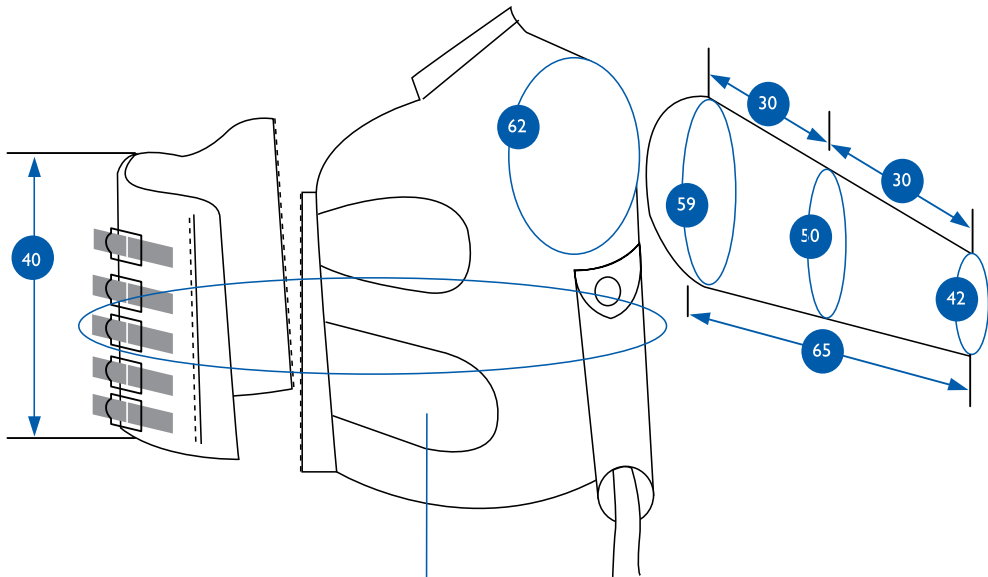

Fastener Sizing			ComfySleeve™		
ITEM	MIN	MAX	SIZE	MIN	MAX
Small Fastener	31	41	S	105	115
Medium Fastener	36	51	M	110	125
Large Fastener	48	66	L	122	140

Arm/Torso Circumference
ONE SIZE
74



12-Chamber Arm Garment

Size Information



12-chamber garments

Overview

Our special design makes arm and shoulder treatment comfortable and easy.

Features


- The 12-chamber arm-shoulder garment provides therapy to the entire shoulder area, as well as the arm.
- Its proprietary design treats right or left arms with a single garment.



Ordering Information

Part No.	Model/Size	Circumference					Length		
		A (Upper Arm)	B	C (Distal End)	Shoulder	Strap (S)	L1	L2	½ L2
L41111C2	LP12 ARM M	63	49	41	76	92-156	91	72	36
L41111C3	LP12 ARM L	74	56	49	88	92-165	106	91	40

* This garment has two zippers, for adjustment of the circumferential fit. The measurements indicated are the maximum (outmost zipper)

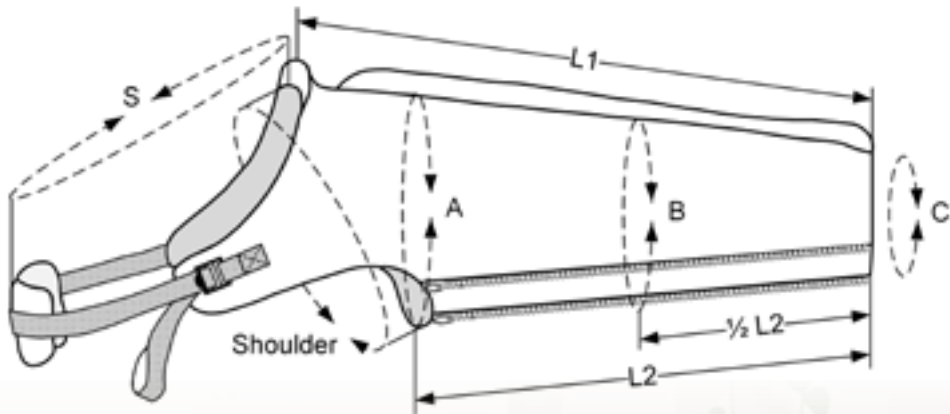

Compatible with



Lympha Press® Optimal Plus



Lympha Press® Mini Press



4/8-Chamber Compression Therapy Systems



Lympha Press® PCD-51™

Overview

Lympha Press® PCD-51™, model 731A offers fully pressure-controlled gradient sequential compression. The compact, user-friendly system features adjustable pressure and time and is designed for quality home treatment.



Features:

- Easy to use and encourages patient compliance
- Compact and easy to carry
- Multi-voltage
- Compatible with Lympha Press® range of unique and innovative garments

Technical specifications:

Cycle	Sequential
Pressure Range	20-80 mmHg
Size	22 cm x 17 cm x 13 cm
Weight	1.5 kg
Electrical	100-240 VAC, 50-60 Hz

Ordering Information

Catalog	Product Name
L10015G3	Lympha Press PCD-51 731A

For indications and contraindications, see page 48



Overview

PCD Pants Lite™ offer an optimal effective treatment of lymphedema in the lower trunk. Its overlapping-chamber compression structure ensures a fluent, comfortable massage and an even pressure application, without gaps.



Features

- Treats the entire lower body, including the abdomen, lower back, hips, and the genital area.
- Adjustable boot accommodates most leg lengths.
- Easy to use – like wearing a pair of pants!
- A zipper to accommodate different sizes.
- Expanders are available


Ordering Information

Model	Cat. No.	Expander	CAT. No
PCD Pants Lite™	L4MPH003	Pants Expanders (pair)	L4Z5LP00
PCD Pants Lite™ Size S	L4MPH005		

Compatible with



PCD-51™

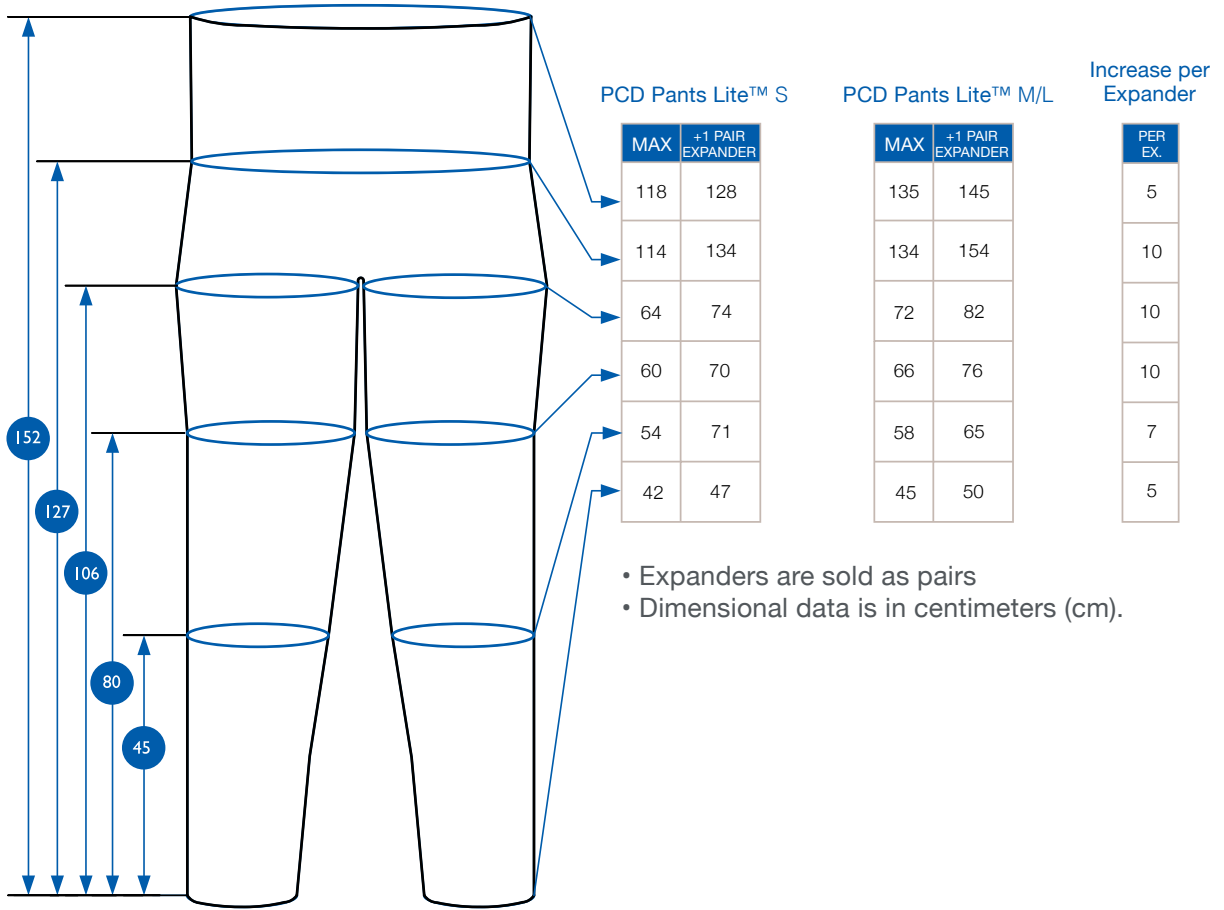

Technical notes

- PCD Pants Lite™ has a single-zipper fastener on each side.



8-Chamber Leg Garment




8-chamber garments

Overview

Ideal for patients who need a garment that is easy to put on and remove.

Features

- Overlapping chambers – our patented system for smooth and effective therapy without backflow or gaps.
- Designed for ease of use.
- Wide selection of standard sizes.
- Specially designed chambers create a conforming fit with air.
- Amputee and custom garments are available.

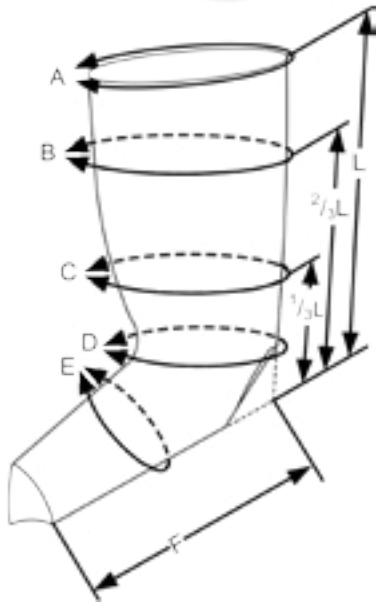

Ordering Information

			Dimensions							With Expander					
Part No	Model/Size	Width	A	B	C	D	E	F	Length	Model	A	B	C	D	E
L41H2365	LP8 3-65	Medium	78	72	60	52	50	35	65	-	-	-	-	-	-
L41H2375	LP8 3-75		80	76	60	52	50	35	75	-	-	-	-	-	-
L41H2385	LP8 3-85		86	80	70	52	50	35	85	-	-	-	-	-	-
L41H2475	LP8 4-75	Wide	102	100	84	66	58	43	75	FT	116	114	97	71	62

Compatible with



PCD-51™

4-Chamber Leg Garment




4-chamber garments

Overview

Ideal for patients who need a garment that is easy to put on and remove.

Features

- Overlapping chambers – our patented system for smooth and effective therapy without backflow or gaps.
- Designed for ease of use.
- Wide selection of standard sizes.
- Specially designed chambers create a conforming fit with air.
- Amputee and custom garments are available.

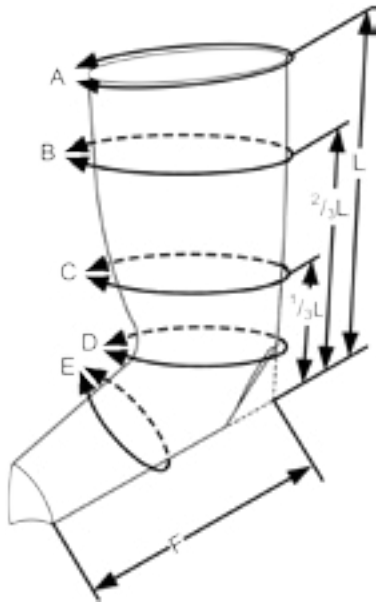

Ordering Information

			Dimensions							With Expander					
Part No	Model/Size	Width	A	B	C	D	E	F	Length	Model	A	B	C	D	E
L41F2250H	LP4 2-50	Narrow	61	54	46	43	48	26	50	AT	71	64	56	50	52
L41F2275H	LP4 2-75		66	56	46	43	48	26	73	-	-	-	-	-	-
L41F2285H	LP4 2-85		70	60	47	43	48	26	84	-	-	-	-	-	-
L41F2295H	LP4 2-95		80	72	50	43	48	26	95	-	-	-	-	-	-
L41F2365H	LP4 3-65	Medium	78	71	65	50	54	38	65	-	-	-	-	-	-
L41F2375H	LP4 3-75		81	74	66	50	50	38	75	-	-	-	-	-	-
L41F2385H	LP4 3-85		86	80	68	53	51	38	85	-	-	-	-	-	-
L41F2465H	LP4 4-65	Wide	105	98	86	67	60	41	65	FT	119	112	99	72	64
L41F2475H	LP4 4-75		107	98	86	72	60	41	75	FT	121	113	99	77	64
L41F2485H	LP4 4-85		109	105	93	67	60	43	85	GT	123	118	106	71	64

Compatible with



PCD-51™


4-Chamber Arm Garment

Overview

Our special design makes arm treatment comfortable and easy.

Features

- Includes an axilla cutout so that the shoulder is treated along with the arm.

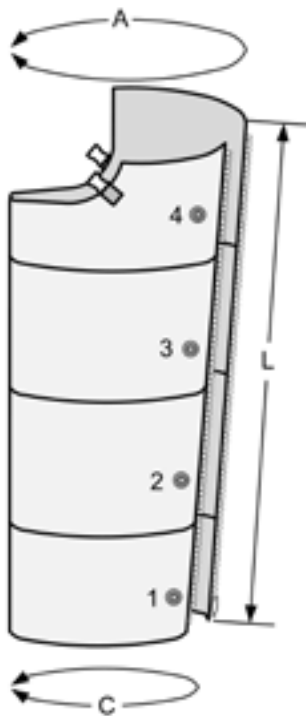

Ordering Information

Part No.	Model/Size	A (Upper Arm)	C (Distal End)	Length
L41F1290H	LP4ARM M/L	58	29	88
L415A087	Expander C	Adds 10 cm to the circumference of the garment		

Compatible with



PCD-51™




Overview

ComfyLite™ treats the entire arm, shoulder, and one side of the back, chest, and abdomen. It provides gentle, targeted therapy to the axilla (armpit), to assist lymphatic drainage.

Features

- Treats arm, chest, shoulder, back, and abdomen.
- Conforms to the patient's body shape, with targeted compression to the armpit area.
- Easy to put on and remove.
- Easy to clean for good hygiene.
- Modular system is cost-effective and easy to use.
- Change the size and side of the garment by switching fasteners.

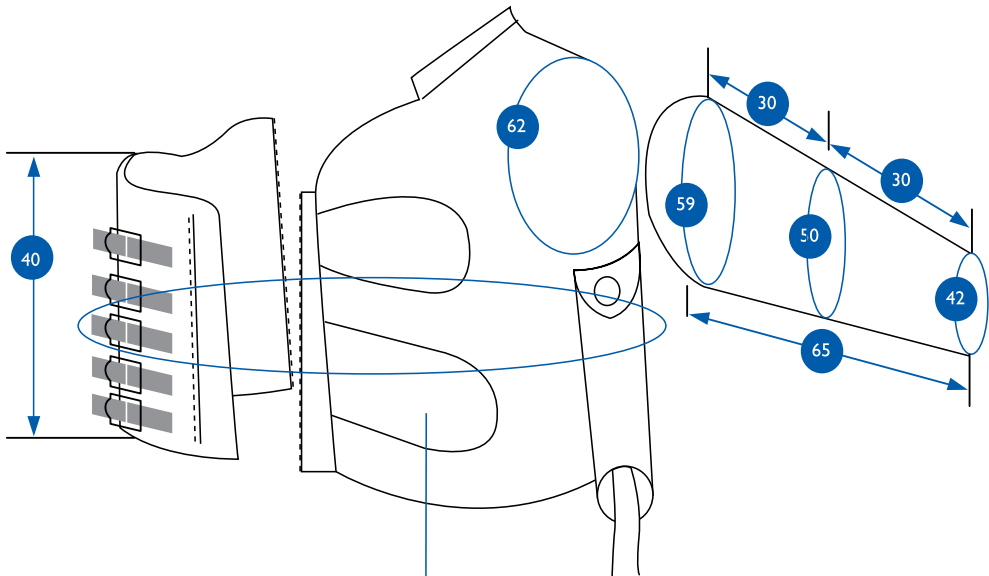

Fastener Sizing		
ITEM	MIN	MAX
Small Fastener	31	41
Medium Fastener	36	51
Large Fastener	48	66

ComfySleeve™		
SIZE	MIN	MAX
S	105	115
M	110	125
L	122	140

Arm/Torso
Circumference

ONE SIZE

74

Ordering Information

MODEL ComfySleeve™ (8 chambers) Use for either arm (not sided)	CAT. No. L41LT100 Arm +Vest
Select: Small Fastener	For Left Arm: L42L0FSL For Right Arm: L42L0FSR
Select: Medium Fastener	For Left Arm: L42L0FML For Right Arm: L42L0FMR
Select: Large Fastener	For Left Arm: L42L0FLL For Right Arm: L42L0FLR

Compatible with



PCD-51™



Garment Fitting Tips

Follow these simple tips to ensure your garment has the best fit for the most effective treatment.



- Measure the patient while the patient is treatment position.
- Obtain all the measurements requested on the Measuring Form. Indicate the position of any unusual protrusions, lobules or landmarks.
- The garments will shorten slightly as they inflate. You may wish to select garment that is slightly longer than the patient's length measurements.
- The garment should not fit tightly at any point. If the patient's circumference is near maximum at any point, choose a wider garment or add an expander insert.

Pants and LymphaPod® Measuring Form

Patient Name _____

Physician Name _____

The diagram illustrates the measurement points for a female figure. It includes a central figure with various measurement lines and points labeled A through G. To the left, there are input boxes for LENGTH, LEFT, and RIGHT measurements. To the right, there are input boxes for CIRCUMFERENCE measurements. At the bottom, there are input boxes for MEASUREMENTS IN: CM and INCHES.

LENGTH

A

B

LEFT RIGHT

C

D

E

F

G

Length Measured from Zero Point

CIRCUMFERENCE

UPPER WAIST

HIPS

GROIN

MID THIGH

KNEE

CALF

ANKLE

INSTEP

FOOT LENGTH

MEASUREMENTS IN:

CM INCHES

Leg Garment Measuring Form

Patient Name _____

Physician Name _____

LEFT RIGHT

A	
B	
C	
D	
E	
F	
G	

Length Measured
from Zero Point

LEFT RIGHT

GORIN	
UPPER 1/3 THIGH	
LOWER 1/3 THIGH	
KNEE	
CALF	
DISTAL GASTROC	
ABOVE ANKLE BONE	
INSTEP	
JOINT LINE	
FOOT LENGTH	

MEASUREMENTS IN:

CM


INCHES

Comfysleeve™ 1-75/ComfyLite™/Lympha Jacket™ Measuring Form

Date _____

Patient Name _____

Physician Name _____

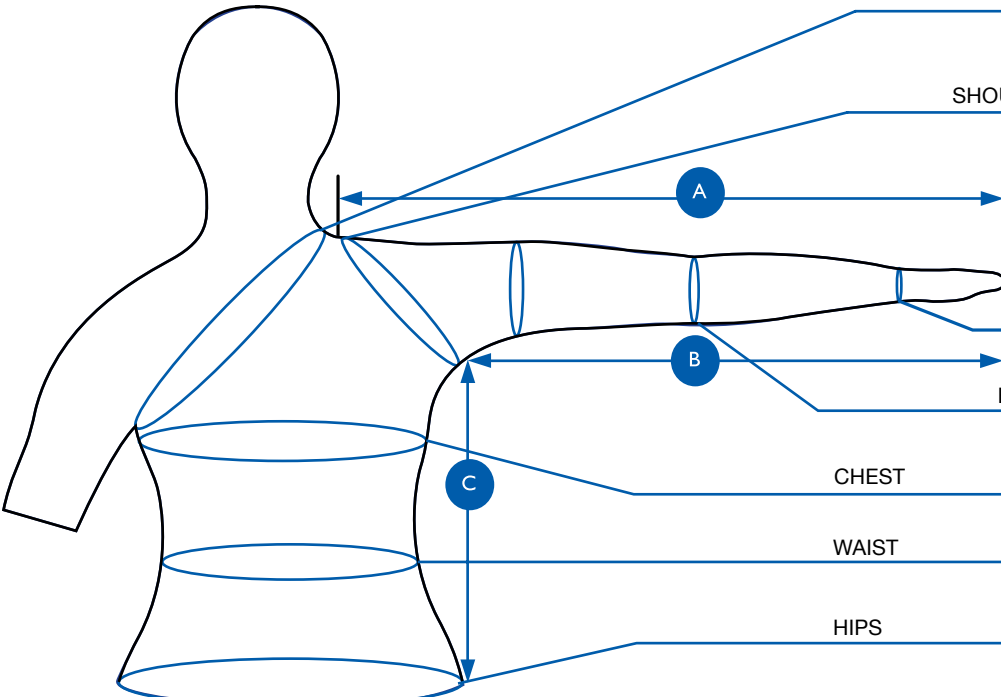

CIRCUMFERENCE

LEFT RIGHT

SHOULDER

SHOULDER

CIRCUMFERENCE

LEFT RIGHT

WRIST

ELBOW

CHEST

WAIST

HIPS

LENGTH A

LEFT RIGHT

LENGTH B

LEFT RIGHT

LENGTH C

LEFT RIGHT

MEASUREMENTS IN:

CM INCHES

Notes _____

44 _____

Arm Garment Measuring Form

Date _____

Patient Name _____

Physician Name _____

LEFT RIGHT

A

B


C

D

E

F

Length Measured from Zero Point

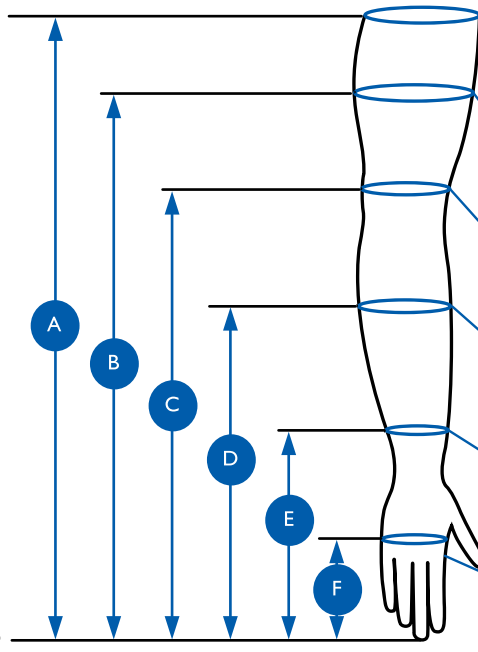

LEFT RIGHT

AXILLA

BICEP

ELBOW

MID FOREARM

WRIST

PALM

Zero Point

MEASUREMENTS IN:

CM INCHES

Notes _____

45 _____

Indications and contraindications

Lympha Press® Optimal Plus

Indications for Use:

Primary lymphedema, secondary lymphedema, venous insufficiency, venous stasis ulcers, dysfunction of the muscle pump, post mastectomy edema, edema following trauma or sports injuries, post immobilization edema. The device may also be beneficial in the management of lipedema, reduction of wound healing time, reduction of pain and swelling after injury and surgery. Lympha Press® Optimal Plus is intended for use by patients at home, and by medical practitioners in clinics.

Contraindications for Use:

Use of Lympha Press® is not recommended in the presence of one or more of the following conditions: known or suspected deep vein thrombosis (DVT) or pulmonary embolism, during the inflammatory phlebitis process, acute infection of the affected limb, decompensated cardiac failure, severe arteriosclerosis or other ischemic vascular disease, any circumstance where increased venous and lymphatic return is undesirable. Due to movement of fluids in the body when using the system, exercise caution when using on patients with heart disease. High pressure is not recommended for patients who have peripheral occlusion disease. The abdominal area should not be treated during pregnancy.

Caution: this device is for sale and use by, or on the order of, a licensed medical practitioner. See the User Manual for additional notes and precautions.

Lympha Press® Mini Press

Indications for Use:

Primary Lymphedema, congenital lymphedema (Milroy's disease), proximal obliteration, distal obliteration, secondary lymphedema, post mastectomy, groin dissection, chronic edema of venous origin, post traumatic edema, lymphatic filariasis. Other indications: lipedema, venous insufficiency, varicose veins, leg ulcers, promotion of wound recovery, postoperative venous ligation or stripping, sport injuries. Intended for use by patients at home, and by medical practitioners in clinics.

Contraindications for Use:

Use of Lympha Press® Mini Press is not recommended in the presence of one or more of the following conditions: sudden pain or numbness, severe arteriosclerosis or other ischemic vascular diseases, unbalanced cardiac insufficiency, known or suspected Deep Vein Thrombosis (DVT) or Thrombophlebitis, gangrene, dermatitis, severe Inflammation of the treated limb, recent skin grafts, untreated or infected wounds. Due to movement of fluids in the body when using the system, exercise caution when using on patients with heart disease. High pressure is not recommended for patients who have peripheral occlusion disease.

The abdominal area should not be treated during pregnancy. **Caution:** this device is for sale and use by, or on the order of, a licensed medical practitioner.

Lympha Press® PCD-51™

Indications for Use:

Primary lymphedema, secondary lymphedema, venous insufficiency, venous stasis ulcers, dysfunction of the muscle pump. PCD-51™ is intended for use by patients at home, and by medical practitioners in clinics.

Contraindications for Use:

Use of the PCD-51™ system is not recommended in the presence of one or more of the following conditions: known or suspected deep vein thrombosis or pulmonary embolism, during the inflammatory phlebitis process, acute infection of the affected limb, decompensated cardiac failure, severe arteriosclerosis or other ischemic vascular disease, any circumstance where increased venous and lymphatic return is undesirable. Due to movement of fluids in the body when using the system, exercise caution when using on patients with heart disease. High pressure is not recommended for patients who have peripheral occlusion disease. **Caution:** this device is for sale and use by, or on the order of, a licensed healthcare practitioner.



Quick Reference Guide

Garments and accessories for use with Lympha Press® Optimal Plus and Lympha Press® Mini Press 12-Chamber Systems

12-Chamber Extremity Garments	See Page 12-29
Comfysleeve® Zipper Leg 2-75	L41L2275
Comfysleeve® Zipper Leg 2-85	L41L2285
Comfysleeve® Zipper Leg 2-95	L41L2295
Comfysleeve® Zipper Leg 3-65	L41L2365
Comfysleeve® Zipper Leg 3-75	L41L2375
Comfysleeve® Zipper Leg 3-85	L41L2385
Comfysleeve® Zipper Leg 4-75	L41L2475
BT Expander (for Comfysleeve 2-75, 3-65)	L415T078
CT Expander (for Comfysleeve 2-85, 3-75)	L415T087
DT Expander (for Comfysleeve 3-85)	.L415T097
ET Expander (for Comfysleeve 2-95)	L415T101
FT Expander (for Comfysleeve 4-75)	L415T465
Lympha Press Arm Size M	L41111C2
Lympha Press Arm Size L	L41111C3
FL Adjustable Leg 2-75	L4132075
FL Adjustable Leg 2-85	L4132085
FL Adjustable Leg 2-95	L4132095
FL Adjustable Leg 2-105	L4132105
FL Adjustable Leg 3-75	L4133075
FL Adjustable Leg 3-85	L4133085
FL Adjustable Leg 3-95	L4133095
FL Adjustable Leg 3-105	L4133105
FL Adjustable Leg 4-75	L4134075
FL Adjustable Leg 4-85	L4134085
FL Adjustable Leg 4-95	L4134095
FL Expander 30X30 cm.	L4Z5A030

Lympha Pants™ and Lympha Pants XL™	See Page 18-19
Lympha Pants™ Hip Hose Exit	L4MP2136C
Lympha Pants™ Heel Hose Exit	L4MP2130C
Lympha Pants™ XL Hip Hose Exit	L4MP2236C
Lympha Pants™ XL Heel Hose Exit	L4MP2230C
Pants Expanders (pair)	L4Z5LP00
LymphaPod®	See Page 20-21
LymphaPod® Left Hose Exit	L4MP800L
LymphaPod® Right Hose Exit	L4MP800R
LymphaPod® Expander Left	L41580AL
LymphaPod® Expander Right	L41580AR
Lympha Jacket™	See Page 24-25
Lympha Jacket™ Short Arm	L4MJ2125
Lympha Jacket™ Long Arm	L4MJ2135
Expander Set for Short Arm Jacket	L4Z52121
Expander Set for Long Arm Jacket	L4Z52130
Jacket Expander for Back	L4150542
Jacket Expander for Arm of Short Arm Jacket	L4150822
Jacket Expander for Arm of Long Arm Jacket	L4150823
Jacket Expander for Back Triangular	L415T510
Jacket Expander for Chest	L4150510
Comfysleeve™ 1-75 Half Jacket Modular System	See Page 26-27
Comfysleeve® 1-75 Left Arm +Vest (no fastener)	L41L0175L
Comfysleeve® 1-75 Right Arm +Vest (no fastener)	L41L0175R
Size S Fastener for Left Arm Treatment	L42L0FSL
Size S Fastener for Right Arm Treatment	L42L0FSR
Size M Fastener for Left Arm Treatment	L42L0FML
Size M Fastener for Right Arm Treatment	L42L0FMR
Size L Fastener for Left Arm Treatment	L42L0FLL
Size L Fastener for Right Arm Treatment	L42L0FLR

Quick Reference Guide

Garments and accessories for use with Lympha Press PCD-51™ 4/8-Chamber Systems

PCD Pants Lite™	See Page 34-35
PCD Pants Lite™	L4MPH003
PCD Pants Lite™ Size S	L4MPH005
PCD Pants™ Expanders (pair)	L4Z5LP00
Leg Garments with 8 Chambers	See Page 36-37
Lympha Press® 8 Leg Size 3-65	L41H2365
Lympha Press® 8 Leg Size 3-75	L41H2375
Lympha Press® 8 Leg Size 3-85	L41H2385
Lympha Press® 8 Leg Size 4-75	L41H2475
Comfylite™ 8 Chamber Half Jacket	See Page 40-41
Comfylite™ Arm +Vest (no fastener)	L41LT100
Size S Fastener for Left Arm Treatment	L42L0FSL
Size S Fastener for Right Arm Treatment	L42L0FSR
Size M Fastener for Left Arm Treatment	L42L0FML
Size M Fastener for Right Arm Treatment	L42L0FMR
Size L Fastener for Left Arm Treatment	L42L0FLL
Size L Fastener for Right Arm Treatment	L42L0FLR

Garments and accessories for use with Lympha Press® PCD-51™ 4/8-Chamber Systems

4 Chamber Arm and Leg Garments	See Pages 30 -41
Lympha Press® 4 Leg Size 2-50	L41F2250H
Lympha Press® 4 Leg Size 2-75	L41F2275H
Lympha Press® 4 Leg Size 2-85	L41F2285H
Lympha Press® 4 Leg Size 2-95	L41F2295H
Lympha Press® 4 Leg Size 3-65	L41F2365H
Lympha Press® 4 Leg Size 3-75	L41F2375H
Lympha Press® 4 Leg Size 3-85	L41F2385H
Lympha Press® 4 Leg Size 4-65	L41F2465H
Lympha Press® 4 Leg Size 4-75	L41F2475H
Lympha Press® 4 Leg Size 4-85	L41F2485H
Lympha Press® 4 Arm Size M/L	L41F1290H
Phlebo Press® Leg Size 2-50	L41P2250H
Phlebo Press® Leg Size 2-75	L41P2275H
Phlebo Press® Leg Size 2-85	L41P2285H
Phlebo Press® Leg Size 2-95	L41P2295H
Phlebo Press® Leg Size 3-65	L41P2365H
Phlebo Press® Leg Size 3-75	L41P2375H
Phlebo Press® Leg Size 3-85	L41P2385H
Phlebo Press® Leg Size 4-65	L41P2465H
Phlebo Press® Leg Size 4-75	L41P2475H
Phlebo Press® Leg Size 4-85	L41P2485H
Phlebo Press® Arm Size M/L	L41P1290H
AT Expander (for 2-50 Leg)	L415T056
FT Expander (for 4-65 and 4-75 Leg)	L415T465
GT Expander (for 4-85 Leg)	L415T485
C Expander (for M/L Arm)	L415A087
FT Expander (for 4-65 and 4-75 Leg)	L415T465

Discover Lympha Press®
World-famous Dynamic
Compression Therapy Systems





Feel the difference. See the results.



info@lymphapress.com | www.lymphapress.com

L991M007