ESTELUX
LUXURY OF AESTHETICS

ESTELUX
LUXURY OF AESTHETICS

Rejuvenel

GLUSKIN

magicalift

WUNDERFILL

Estelux
Filler

BIOHA
LABORATORIES

ESTELUX ŞİRKƏTİ HAQQINDA

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BIOHA
LABORATORIES

ABOUT ESTELUX COMPANY

Estelux company has been officially operating in the Republic of Azerbaijan since 2018. Our company is currently the exclusive distributor of 6 international brands and 59 product types and collaborates with 1700 cosmetologists. We have stands in the international arena at Amwc Monaco, Dubai Derma, Imcas Paris, Imcas Shanghai, Kolkhida Tbilisi, Beauty Istanbul, Cosmoderma, and at local exhibitions in our country such as Baku Derma and others. The products we distribute have ISO, CE, GMP certificates and licenses from the Ministry of Health. We import cosmetology and dermatology products to our country from United Kingdom, Turkey, Netherlands, Korea, and China.

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ESTELUX
LUXURY OF AESTHETICS

PROCEDURE INSTRUCTIONS

Before the treatment, makeup and all other cosmetic products must be removed to ensure the skin is clean. This can be remedied by using a facial cleansing lotion. Using cold water will ensure that pores do not expand. Medications such as Aspirin, Advil, etc. should be taken at least 3 days before the treatment. Because they can increase swelling or bleeding at the injection site. It is also recommended to stay away from alcohol, fish oil, ginkgo, and omegas before the treatment. During the first day after the cosmetology and dermatology treatment, try to keep the face relaxed by not making any highly expressive facial movements and avoiding any strenuous exercise. It is also important not to be exposed to excessive heat (e.g., sauna, solarium, etc.). Because this can reduce the longevity of the results. Do not massage the injected area unless instructed. Makeup can be reapplied immediately after the treatment. It is important to remember to avoid the following things before and after treatment: Advil, aspirin or other anti-inflammatory drugs, Alcohol, Fish oil, Omegas, Strenuous exercise, Excessive heat or cold.

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POLYLACTIC ACID FILLER / POLİLAKTİK DOLGU

Product Descriptions

Function: A biostimulator for injection made of Poly-L-lactic acid. When injected into the deep dermis or subcutaneous tissue, PLLA creates an immediate filling effect. Over time, PLLA begins to break down. It stimulates skin tissue fibroblasts, increases collagen, and has an elastic plumping effect.

Composition: Poly-L-lactic acid microspheres (PLLA) 150 mg, Cross-linked hyaluronic acid (17mg/mL), Lidocaine hydrochloride (3mg/mL), Phosphate buffer.

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
CALCIUM FILLER / KALSİYUM DOLGU

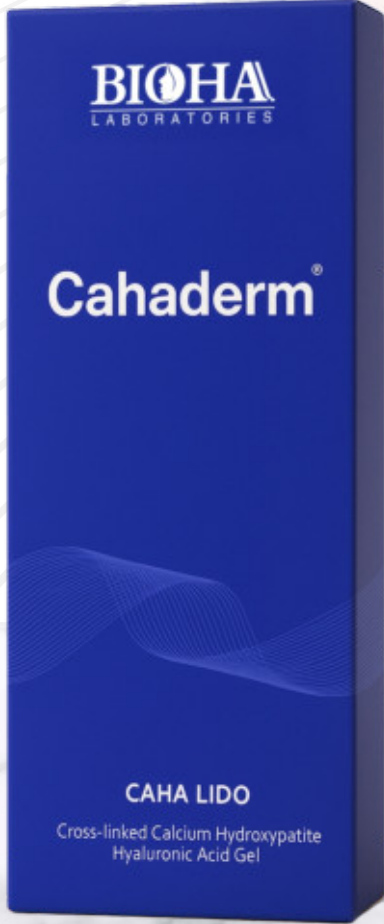Product Descriptions

Function: Stimulates the restoration of collagen and elastin. Restores the correction of facial soft tissue deficiency. It is placed under the skin for deep and superficial biostimulation. It is injected with a cannula in the jawline and pre-auricular area.

Composition: 30% Calcium hydroxyapatite microspheres, 25-45 micron diameter, 70% Cross-linked sodium hyaluronate gel (20mg/ml), Phosphate buffer, Lidocaine hydrochloride (3mg/ml).

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5-POINT LIFTING / 5-NOKTA LİFTİNG

Product Descriptions

Function: Strengthens, moisturizes, and brightens the skin. It restores facial sagging and neck lines, is anti-aging, and gives a 5-point lifting effect. It balances interstitial moisture, evens out skin tone, restores the sensitive skin barrier, reduces wrinkles, and lessens the effects of discomforts like eczema and psoriasis.

Composition: 32mg(H-HA) + 32mg(L-HA) / 2ml.

Protocol: Inject from 5 points on each side of the face, 0.2 ml per point. On the neck, 10 points, 0.2 ml per point ; a total of 4 ml for the face and neck.

Hyaluronic Acid Density: 32 mg / ml.

Needle Size: 27G, 30G.

Duration: 6-18 months.

Storage Temperature: 2-30°C.

Medical cross-linked sodium hyaluronate gel soft tissue filler.

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BIOREVITALIZATION / BİOREVİTALİZASYON

Product Descriptions

Function: Moisturizes and nourishes especially sensitive, dry skin and helps in solving the problem of spots. It is recommended for skin repair after sunburn. Macromolecules retain moisture for a long time, small molecules moisturize quickly, and micro molecules penetrate the dermis to moisturize and bind water. It eliminates free radicals and promotes cell redox, preserves the vitality of the skin, and helps to reverse the aging process.

Composition: High Molecular HA, Medium Molecule HA-TUM/HA-T, Acetylated HA, Low Molecular HA-TLM, MINI Hydrolyzed HA.

Protocol: 1-2 vials for one treatment, injected into the upper dermis or at the junction of the epidermis and dermis. A 2-3 week injection interval is recommended. A full course of treatment is five times.

Phase One: Three intensive treatments (The skin becomes moisturized and radiant with tightened pores.).

Phase Two: Twice with a 1-month interval (Fibroblasts are restored and the skin becomes elastic.).

Phase Three: Once every two months to enhance the effect (Fibroblasts continue to secrete collagen and elastin to make the skin shiny and firm, reducing dry and fine lines.).

Can be used with Meso Needle, Meso Gun, Derma Roller, and Derma Pen.

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Reyoungel®
MESO LIFTING
SKIN REJUVENATION SOLUTION

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LIPOLYTIC / LİPOLİTİK

Product Descriptions

Function: Helps in dissolving excess fats. The solution reduces the formation of cellulite and promotes the burning of fats. It can be injected into fat tissue. It can be injected into the face, buttocks, double chin, neck, chin, abdomen, waist, back, arms, hands, legs. Contains DMAE.

Composition: Hyaluronic acid, Nicotinamide, Deoxycholic acid, DMAE 5 mg.

Protocol: Cheeks, chin, neck: one-time treatment 3-4ml, 0.2-0.4ml for each injection site. The distance of the injection site cannot be less than 1 cm, and the depth cannot be more than 0.5 cm. A 2-week injection interval is recommended. Abdomen, waist, arms, legs, hips: one-time treatment 5ml, 0.5-0.8 ml for each injection site. The distance of the injection site cannot be less than 1 cm, and the depth cannot be more than 1 cm. A 1-week injection interval is recommended. You can see the result of the treatment after 4-6 times. It is necessary to monitor the treatment according to the patients.

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PARIS GLOW / PARİS İŞİLTİSİ

Product Descriptions

Function: It is a "Paris glow" mesotherapy. It moisturizes the skin, tightens pores, smooths, and supports collagen. It creates a magnificent radiance and vitality. It strengthens the skin's elasticity, and ensures the restoration of dry and sensitive skin that has lost its balance and radiance. It reduces skin stress and increases its quality. It is an anti-inflammatory. It accelerates cell metabolism.

Composition: Hyaluronic Acid, Mannitol, Arginine Acetyl Hexapeptide-8, Panthenol.

Protocol: To biologically rejuvenate the skin: Micropomfus technique, 0.05-0.1 ml for each pomfus. The number of sessions is determined every 21 days depending on the skin's reaction to the treatment. For maintenance, 1 session every 6 months is recommended. The maximum recommended dose in a single session is 2.5 ml. It is recommended to use a 1ml locked syringe attached to a 30G 13mm – 4mm needle.

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6-POINT LIFTING / 6-NOKTA LİFTİNG

Product Descriptions

Function: It is a six-point lifting miracle. It forms the skin's moisture balance and gives a magnificent lifting effect. It increases the effectiveness of the lifting by renewing elastin and collagen fibers. The skin tone is lightened and evened out over the entire face.

Composition: Hyaluronic Acid, Succinic Acid, Glutathione, Arginine, Cystine, Glycine, Glutamic Acid, Tryptophan, Oligopeptide, Aspartic Acid, Pyridoxine, Cyanocobalamin, Ascorbic Acid, Tocopheryl Acetate, Panthenol, Thiaminbiotin.

Protocol: The biological lifting solution is injected into the skin with a 6-point technique. It is carried out by injecting into 6 different points determined specifically for the patient on both sides of the face, for a total of 12. 0.2 ml per point. A lifting effect is seen in a short time. The number of sessions is determined every 21 days depending on the skin's reaction to the treatment. For maintenance, 1 session every 6 months is recommended. The maximum recommended dose in a single session is 2.4 ml. It is recommended to use a 1ml locked syringe attached to a 30G 13mm – 4mm needle.

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Product Descriptions

Function: Removes periorbital and under-eye wrinkles. It helps to eliminate bruises, dark circles, and other unpleasant appearances. It consists of 8 different B complex vitamins. It has strong antioxidant properties and increases red blood cells. It contains Hexapeptide to remove fine wrinkles.

Composition: Hyaluronic Acid, Thiamine, Ascorbic Acid, B-complex, Succinic Acid, Vitamin K, Cobalamin, Mannitol, Panthenol, L-arginine, Acetyl Hexapeptide-8, Glutathione, Riboflavin.

Protocol: It is injected with the Micropomfus technique at a depth of 1-1.5 mm under the eye and 1.5-2 mm around the eye. 0.1 ml for each pomfus. It should be applied in sessions every 15-21 days. For maintenance, 1 session every 6 months is recommended. The maximum recommended dose in a single session is 1.5 ml. (for a total of 2 eyes and their surroundings) 0.3 ml under each eye and 0.2 ml around it. It is recommended to use a 1 ml locked syringe attached to a 30 G 4 mm needle.

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HAIR MESOTHERAPY / SAÇ MEZOSU

Product Descriptions

Function: Prevents hair loss, lengthens small hairs, supports vitamins at the hair roots, softens the hair, and gives it vitality. It strengthens the hair roots. Since it contains biotin and panthenol, it prevents oiliness.

Composition: Hyaluronic Acid, Biotin, Ascorbic Acid, Niacinamide, Panthenol, Tocopheryl Acetate, Pyridoxine.

Protocol: It is injected with the Micropomfus technique to biologically rejuvenate the scalp. 0.1 ml for each pomfus. The number of sessions is determined every 10-15 days depending on the scalp's reaction to the treatment. Results can be achieved in 3 sessions. For maintenance, 1 session every 6 months is recommended. The maximum recommended dose for the entire head area in a single session is 5 ml. It is recommended to use a 1 ml locked syringe attached to a 30 G 4 mm needle.

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SPOT MESOTHERAPY / LEKE MEZOSU

Product Descriptions

Function: Ensures the removal of sun, laser, pregnancy, childbirth, and other specific spots. It restores the skin after acne. It ensures the pigmentation in the skin. It prevents the formation of spots by reducing the tyrosinase enzyme.

Composition: Hyaluronic Acid, Panthenol, Ascorbic Acid, Glutathione, Niacinamide, Riboflavin, Thiamine, Cyanocobalamin, Retinol, Succinic Acid, Arbutin, Tranexamic Acid, Mannitol, Rutin.

Protocol: It is injected with the Micropomfus technique into areas with different skin color and all over the face. 0.1 ml for each pomfus. The number of sessions is determined every 15 days depending on the skin's reaction to the treatment. The maximum dose to be applied in one session is 2.5 ml. It is recommended to use a 1ml locked syringe attached to a 30G 13mm – 4mm needle.

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MOISTURIZING TREATMENT / NEMLENDİRİCİ

Product Descriptions

Function: Hyaluronic acid attracts moisture to the skin like a magnet and retains it. It supports and strengthens the skin's natural protective barrier. It softens, smooths, moisturizes dry skin, and gives it elasticity and radiance. It smooths wrinkles and increases the production of new collagen. It is the main component of the intercellular matrix.

Composition: Purified Hyaluronic Acid.

Protocol: It is injected with the Micropomfus technique to biologically rejuvenate and moisturize the skin. 0.05-0.1 ml for each pomfus. For filling wrinkles: Linear technique, 0.1 ml – 0.2 ml per wrinkle. The number of sessions is determined every 21 days depending on how the skin reacts to the treatment. For maintenance, 1 session every 6 months is recommended. The maximum recommended dose in a single session is 2.5 ml. It is recommended to use a 1ml locked syringe attached to a 30G 13mm – 4mm needle.

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COLLAGEN REGENERATION / KOLLAJEN
REJENERASYONU

Product Descriptions

Function: Collagen renews the dermis and connective tissues. It increases the skin's elasticity, reduces the signs of aging and wrinkles. It removes scars and acne marks. It restores the elasticity in the skin and enhances vitamins. Collagen is perfectly provided.

Composition: Purified Collagen.

Protocol: It is injected with the Micropomfus technique to biologically rejuvenate the skin and restore collagen. 0.05-0.1 ml for each pomfus. For filling wrinkles: Linear technique, 0.1 ml – 0.2 ml per wrinkle. The number of sessions is determined every 21 days depending on how the skin reacts to the treatment. For maintenance, 1 session every 6 months is recommended. The maximum recommended dose in a single session is 2.5 ml.

It is recommended to use a 1ml locked syringe attached to a 30G 13mm – 4mm needle.

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BLUE SHINE

5ML

STRENGTHENING TREATMENT / GÜÇLENDİRİCİ

Product Descriptions

Function: It is a skin-strengthening mesotherapy. It moisturizes the lower layer of the skin, supports skin health because there are different vitamins, and activates the activity of cells. It nourishes the skin and increases the effectiveness of peptides. Amino acids increase the elasticity of the skin, reduce sagging and wrinkles, and nourish the skin. It tightens pores, supports collagen, and moisturizes. It removes fine wrinkles in the skin, providing a magnificent radiance.

Composition: Hyaluronic Acid, Aspartic Acid, Arginine, Lysine, Niacinamide, Cystine, Glutathione, Glycine, Mannitol, Oligopeptide-3, Tranexamic Acid, L-tryptophan, Sodium Succinate, Glutamic Acid, Salmon Dna, 24 Amino Acids, 14 Vitamins.

Protocol: It is injected with the Micropomfus technique to biologically rejuvenate the skin. 0.1 ml for each pomfus. The number of sessions is determined every 15 days depending on the skin's reaction to the treatment. For maintenance, 1 session every 3 months is recommended. The maximum dose to be applied in one session is 2.5 ml. It is recommended to use a 1ml locked syringe attached to a 30G 13mm – 4mm needle.

Since Gloskin mesotherapies are produced with new generation damping technology, no pain, burning, or stinging is felt during the injection.

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PINK COCKTAIL / PEMBE KOKTEYL

Product Descriptions

Function: Gives the skin instant lifting and radiance. It moisturizes and nourishes the skin. It has a strong anti-aging effect. It is a pink cocktail enriched with a vitamin complex. It plays a role in the production of red blood cells, thus eliminating the pale and yellow appearance of the skin. It protects the skin and prevents the formation of wrinkles by balancing sebum production. It helps to maintain a youthful appearance.

Composition: Hyaluronic Acid, Cyanocobalamin, Niacinamide, Ascorbic Acid, Arginine, Cystine, Glycine, Glutamic Acid, Tryptophan, Glutathione, Vitamins, Salmon Dna.

Protocol: It is injected with the Micropomfus technique to biologically rejuvenate and moisturize the skin. 0.05-0.1 ml for each pomfus. The number of sessions is determined every 15-21 days depending on how the skin reacts to the treatment. For maintenance, 1 session every 6 months is recommended. The maximum recommended dose in a single session is 5 ml. It is recommended to use a 1ml locked syringe attached to a 30G 13mm – 4mm needle.

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ACNE MESOTHERAPY / AKNE MEZOSU

Product Descriptions

Function: Treats acne on the skin. It prevents breakouts, moisturizes, and directly treats acne and inflammatory breakouts on the skin. It creates a brighter and more vibrant skin barrier.

Composition: Hyaluronic Acid, Glycolic Acid, Salicylic Acid, Vitamin A.

Protocol: It is injected with the Micropomfus technique for treating acne on the skin. 0.05-0.1 ml for each pomfus. The number of sessions is determined every 15-21 days depending on how the skin reacts to the treatment. For maintenance, 1 session every 6 months is recommended. The maximum recommended dose in a single session is 2.5 ml. It is recommended to use a 1ml locked syringe attached to a 30G 13mm – 4mm needle.

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9-POINT LIFTING / 9-NOKTA LİFTİNG

Product Descriptions

Function: It is a 9-point lifting miracle. It forms the skin's moisture balance and gives a magnificent lifting effect. It increases the effectiveness of the lifting by renewing elastin and collagen fibers. The skin tone is lightened and evened out over the entire face.

Composition: Hyaluronic acid, Succinic acid, Glutathione, L-arginine, Niacinamide, Glycine, Glutamic acid, Tryptophan, Oligopeptide-1, D-Aspartic acid, Pyridoxine, Cyanocobalamin, Ascorbic acid, Tocopherol acetate, Thiamine HLC, biotin, zinc, Calcium chloride, Potassium chloride, Magnesium chloride.

Protocol: The biological lifting solution is injected into the skin with a 9-point technique. It is carried out by injecting into 9 different points determined specifically for the patient on both sides of the face, for a total of 18. 0.15 ml per point. A lifting effect is seen in a short time. The number of sessions is determined every 21 days depending on the skin's reaction to the treatment. For maintenance, 1 session every 6 months is recommended. The maximum recommended dose in a single session is 2.7 ml. It is recommended to use a 1ml locked syringe attached to a 30G 13mm – 4mm needle.

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ADVANTAGES OF TECHNOLOGY

CROSS-LINKING TECHNOLOGY

“SAX-HA” technology is a new type of cross-linking technology that allows synthetic hydrogels with hydrophobicity in the polymer chain to achieve strong self-association functions through strong hydrophobic association. This self-association is dynamic, and the injected gels tend to stick together and aggregate to maintain shape and reduce gel migration in the body.

There are two commonly used crosslinking agents: BDDE (butylene glycol diglycidyl ether) and DVS (divinyl sulfone). They have certain residues and adverse reaction reports.

crosslinking system
[Low cross-linking efficiency]

HCXL system
[High cross-linking efficiency]
Hyaluronic acid BDDE (crosslinker) Binding site
Non-binding site

Reyoungel uses SAX-HA cross-linking technology, which has the following advantages:

1. More thorough removal of cross-linking agent (BDDE) residues, reducing adverse reactions.
2. Better retention, reduce displacement. Make the shaping effect more natural and long-lasting (6-18 months).

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AB Tasarım İnceleme sertifikası