

AnyreachC™ & AnyreachP™

PTCA Guidewire

ANYREACH, ANYWHERE

APT Medical provides a series of **Anyreach PTCA guidewire** features with a Safe Net Core and Core-to-Tip design to offer superior torque control and safety, making you more confident in your clinical practice.



Support shaft: stainless steel

Proximal stainless steel core with PTFE coating provides excellent pushability and rail support

Core-to-tip design

Core-to-tip to enhance the tactile feedback and pushability



Hydrophilic coating

Durable hydrophilic coating ensures low friction advancement


Safe Net Core Design

Offers superior torque control, safety and shape retention

AnyreachC™

C for Coil: the metal coil gives you real tactile feedback

Metal coil, safe net and Nitinol core combination introduces a unique structure, maximize the torqueability and safety, different tip loads available to make it an optimized choice for your everyday practice.

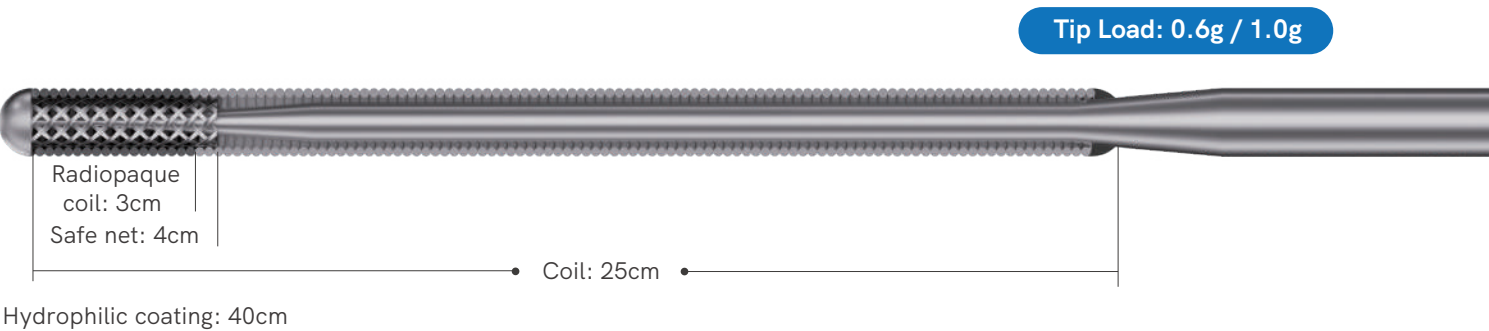

The lightest tip load of 0.3g is available for retrograde approach


Tip Load: 0.3g

Safe net length: 4cm
Stainless steel coil: 16cm
Hydrophilic coating: 55cm


AnyreachP™

P for Polymer: the smooth feel you wish for

Polymer sleeve on the distal covered with hydrophilic coating for enhanced lubricity to provide excellent trackability and crossability in tortuous vessels.

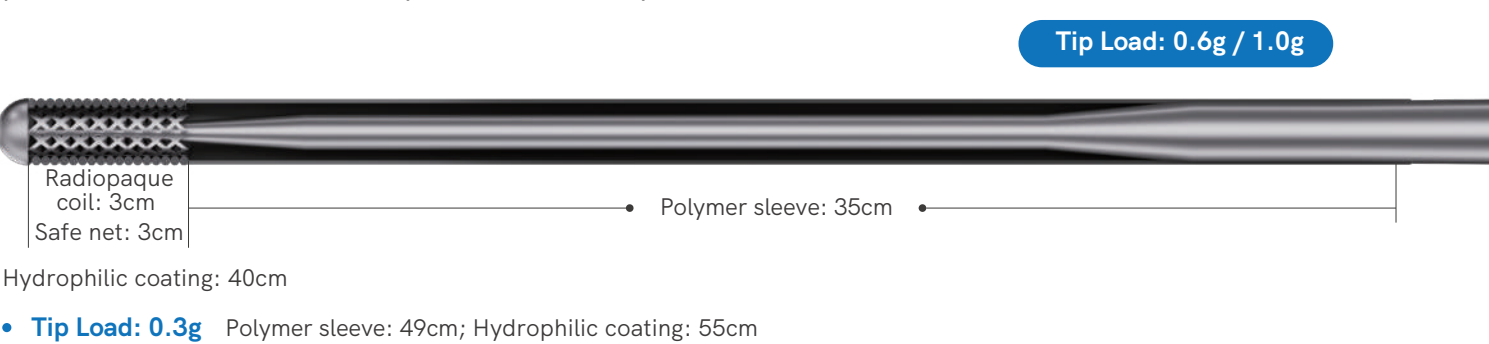


AnyreachP™ Lance℞ (CTO Indicated)

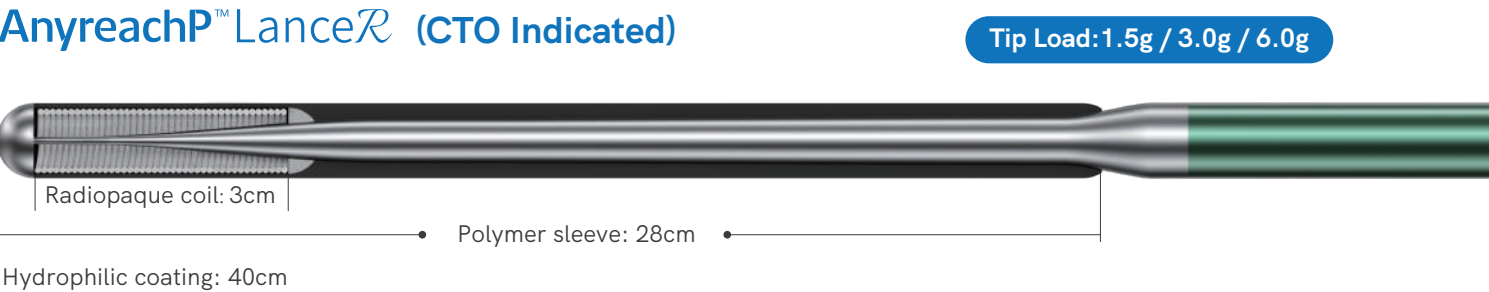

ORDERING INFORMATION

Workhorse

Product	REF	Tip Load	Curve	Length
AnyreachC™	40S019B06	0.6g	Straight	190cm
	40S030B06	0.6g	Straight	300cm
	40A019B0634	0.6g	A (3mm 45°)	190cm
	40H019B0634	0.6g	H (3mm 45°)	190cm
	40S019C06	1.0g	Straight	190cm
	40S030C06	1.0g	Straight	300cm
	40A019C0634	1.0g	A (3mm 45°)	190cm
	40H019C0634	1.0g	H (3mm 45°)	190cm
Product	REF	Tip Load	Curve	Length
AnyreachP™	59S019B0636	0.6g	Straight	190cm
	59S030B0636	0.6g	Straight	300cm
	59A019B063634	0.6g	A (3mm 45°)	190cm
	59H019B063634	0.6g	H (3mm 45°)	190cm
	59S019C0636	1.0g	Straight	190cm
	59S030C0636	1.0g	Straight	300cm
	59A019C063634	1.0g	A (3mm 45°)	190cm
	59H019C063634	1.0g	H (3mm 45°)	190cm

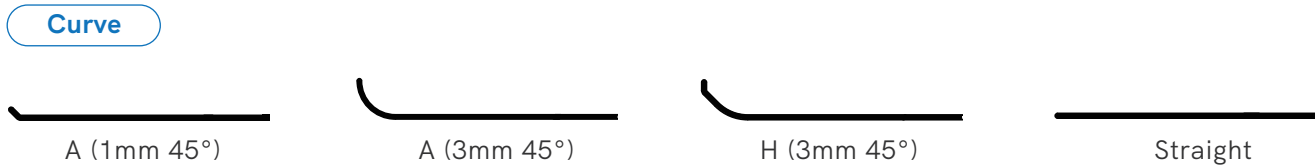

Workhorse + Retrograde

Product	REF	Tip Load	Curve	Length
AnyreachC™	40A019C0614	0.3g	A (1mm 45°)	190cm
	40A030C0614	0.3g	A (1mm 45°)	300cm
AnyreachP™	59A019C063614	0.3g	A (1mm 45°)	190cm
	59A030C063614	0.3g	A (1mm 45°)	300cm

AnyreachP™ Lance℞ (CTO Indicated)

Product	REF	Tip Load	Curve	Length	Color of tubeclip
AnyreachP™ Lance℞	59H119B083634	1.5g	H (3mm 45°)	190cm	●
	59S019B0836	1.5g	Straight	190cm	
	59H119C083634	3.0g	H (3mm 45°)	190cm	●
	59S019C0836	3.0g	Straight	190cm	
	59S119C1036	6.0g	Straight	190cm	●
	59H119C103634	6.0g	H (3mm 45°)	190cm	

For LanceR series, different tip loads will be indicated by different colors of tubeclip.

Indicators for models identification (at proximal end)

