



Beijing ADSS Development Co.,Ltd.

Headquarter address:

F9, Greenland Wealth Plaza, Jinxing Road, Daxing District, Beijing, China

Tel: 0086-10-83625120 Fax: 0086-10-83625121

Email: adss@adss.com.cn

Factory address:

An De Industrial Zone, FuDa Road(East side), TongSheng Road(North side), Gu'an Industrial Park SouthZone, Lang Fang City, Hebei Province, China 102600



Version No.: 2023-12 Designer: Vivian Zhou Verifier: Monica Shao

ADSS®

COMPANY

PROFILE

Beijing ADSS Development Co., Ltd., established in July 2005, is a professional manufacturer of Aesthetic & Medical Laser Device. Now ADSS has over 500 employees and 160,000m2 self-owned industrial park which covers R&D Center, Testing Lab, Automatic Machining Workshop, Production Lines, Service Support, Training School, etc. ADSS can match different market requirements for OEM&ODM, and professional clinical guide.

ADSS always sticks to the principle of "Safety, Effectiveness, Reliability" over the years, which helps ADSS achieve good reputation in more than 140 countries. Being the winner of the "Most Outstanding Chinese Entry" company on Bologna Cosmoprof, ADSS is your one-step cooperation partner!

INDEPENDENT R&D



ADSS has its own R&D center with independent R&D capabilities. The R&D team masters the core technology of all products. Besides, ADSS still insists on R&D investment and keeps the increasing every year so that our technology has always been at the leading position in China.

2 COMPLETE INDUSTRYL CHAIN



ADSS has imported machine tools from World-leading brands, such as US Haas, Japan FANUC, Italy Salvagnini and Prima Power, to reach achievement of Whole Industry Chain and also ensure High-quality Manufacturing, Fast Delivery and Quality Control. ADSS is committed to providing the safest, and the most effective high-quality equipment & services to exceed customer expectations.



3 SUPPORT OEM/ODM SERVICES



ADSS can provide customers with one-stop services from product development, design assembling, packaging, etc to the most, cost-effective solutions.

4 AFTER-SALES SERVICE



From 2005 to now, ADSS has had rich experience in international after-sales service. All machines are designed and improved with optimal internal structure to make maitnenance easy. 24 hours online service support is also available.

5 COMPANY CULTURE



ADSS has its own corporate philosopy and management philosopy, and puts great importance on corporate culture. Culture is faith. Only with common faith can ADSS develop bigger and build long-term relationships with clients.

(E ISO 13485 ROHS NMPA FDA

ADSS HISTORY

2013-2017

Business Expanded

- Moved into self-owned headoffice building
- Started Phase1 construction of ADSS Industrial Park
- Increased R&D investment
- Obtained US FDA certificate
- Obtained CFDA certificate of Class III laser

2021-2022

Leading Aesthetic & Medial Device Manufacturer

- > Started Phrase 1 construction of ADSS Industrial Park
- > Expanded R&D team to 60 engineers
- Imported machine tools from USA and Japan
- Established comprehensive researching and testing center with China TOP 1 Hospital PUMCH













2005-2012

Company Establishment & Growth

- > Established in July 2005
- ➤ Obtained ISO9001
- ➤ Obtained ISO13485
- ➤ Obtained 1st Germany TUV Medical CE

2018-2020

Production Capability Upgraded

- Moved into 160,000m² self-owned Industrial Park
- > Set up largest dust-free workshop for optical laser production
- Built testing labs for EMC and LVD
- ➤ Being Top 1 beauty machine exporter in China

2023-2024

Whole Industry Chain Layout

- > Imported machine tools from Italy Salvagnini and Prima Power
- > 2023 expend CNC Industrial chain and production workshop
- > ADSS has developed higher technology and added more machine (2024)
- Got approval to build the second Industrial Park in Beijing Biomedical Base, which covers an area of 20,000m2 (2024)







PRODUCT

HAIR REMOVAL SERIES	01
Diode Laser Hair Removal Machine - 4D TECLASER 1064nm Long Pulse	(01-04)
INTENSE PULSED LIGHT SERIES	02
Multifunctional DPL	(05-08)
BODY SLIMMING SERIES	03
Ande Smart 360°Cryolipolysis Machine - COOL 360 Muscle Building Machine - EM Contouring / EM MAX 1060NM Diode LaserShape Machine Multishape Body Slimming Machine - VS+/VS++ Legenda Deep Beauty 448KHz	(09-22)
HIFEM PELVIC FLOOR MUSCLES CHAIR	04
ME Bella / ME Bella Pro	(23-24)

PIGMENTED LESIONS REMOVAL SERIES	05
Q-Switched Nd:YAG Laser	
Super Picosecond Laser	(25-28)
	(23 20)
LASER SKIN RESURFACING SERIES	06
UltraPulse Fractional CO2 Laser	
Golden Age Microneedling RF	(20.72)
	(29-32)
ANTI-AGING SERIES	07
ME Time	
360 Ultra HIFU (High Intensity Focused Ultrasound)	
20MHz Cold CRF	
UltraTone Multifunctional 2 IN 1 System	(33-40)
FACIAL SKIN CARE SERIES	08
All In One Hydrafacial	
Bio Light Skin Repairing System	
HydraMagic Mesotherapy Injection System 3D Skin Analyzer - ViewSkin	
	(41-48)

4D TECLASER (2400W / 1600W / 600W) DIODE LASER HAIR REMOVAL SYSTEM

The most advanced diode laser hair removal system that combines 4 wavelengths to permanently and safely remove any type of hair on all skin types.

The high energy density and advanced cooling system make the hair removal procedure much faster and pain free. The device leverages the benefits of 4 combined wavelengths 755nm & 808nm & 940nm & 1064nm, with an improved patient experience and a business-oriented approach, creating a new solution in the world of professional hair removal.











WORKING THEORY


4D TECLASER hair removal System, adopts the light selective absorption theory, heating the hair shaft and hair follicle to destroy the hair follicle and oxygen organization around hair follicle, thereby to get the permanent hair removal purpose.

755nm Light-colored and thin hair.

808nm Golden standard for all skin types and hair colors.

940nm coagulates the microvessels around the hair follicle, cutting off the nutrient supply of hair follicles.

1064nm A perfect solution for darker skin types.

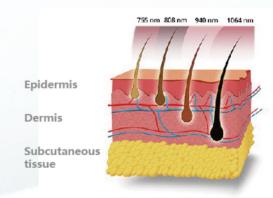

INTUITIVE AND USER-FRIENDLY INTERFACE

Customized for male and female.
For all skin types I - VI.
5 cooling levels adjustable.
Multi languages menu.





HANDPIECE

Handpiece comes with a smart touched screen and parameters can be displayed and set freely on handpiece as well.



FACIAL TIP (Optional)

Easily treat areas that are normally hard to reach, including the ears, nostrils and glabella.

BEFORE & AFTER





SPECIFICATIONS

Frequency	1-10Hz adjustable
Cooling	Air + Water + TEC
Laser Type	755nm/808nm/940nm/1064nm

1064NM LONG PULSE LASER HAIR REMOVAL

The most effective and safest laser for dark skin types.



Function

✓ Hair Removal


▼ Skin Rejuvenation

✓ Vascular Removal



WORKING THEORY

Long Pulsed ND:YAG is based on the principle of photo-thermal destruction of the target tissue. Target structures (blood vessels and pigmented lesions) absorb the focussed Laser light, (which is set at a specific fluence and pulse width based on the condition being treated). Light energy gets converted to thermal energy (heat), thus destroying the target structure.

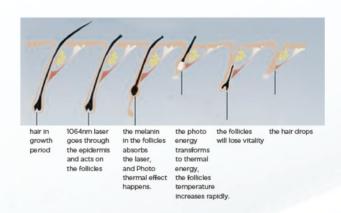

BEFORE & AFTER



HANDPIECE

Japanese Mistsubishi fiber laser, high power and stable output.
Skin contact cooling.





INTERFACE





SPECIFICATIONS

Screen size	10.4"
Wavelength	1064nm
Frequency	1-10Hz adjustable
Spot size	2-10mm adjustable
Dimension of machine	W33*D86*H92cm
Weight	78kg

MULTIFUNCTIONAL DPL

DPL technology a new generation of IPL (Intense Pulsed Light), for multifunctional application like permanent hair removal and skin rejuvenation, which is having a sweeping success.



Function

Mair removal

✓ Skin rejuvenation & whitening

✓ Acne scars removal

Vascular removal

✓ Pigments removal



HISTORY OF INTENSE PULSED LIGHT


IPL (Intense Pulsed Light) describes the use of intense pulses of non-coherent light distributed over a range of wavelengths. The technology utilizes multipulses of energy at high intensity to target various chromophores in the skin.

Wavelength: 400nm - 1200nm

E-light (IPL + RF) RF is a noninvasive technique that raises the temperature of the deep layers of the skin. Thus, RF complements IPL, expanding the effect. Wavelength: 400nm - 1200nm

SHR (Super Hair Removal) uses a single pulse energy and a high repetition rate, which allows the skin to be heated gently and the energy is then accumulated by the skin. With the In-Motion technology (up to 10 shots/second), make the procedure more efficient. Wavelength: 530nm - 950nm; 640nm-950nm


DPL the spectrum of intense pulsed light is compressed in 530nm - 750nm. Its particular narrow-spectrum light can improve the ratio of absorption rate of melanin and hemoglobin, which results in better performance of hair removal and vascular theapy.

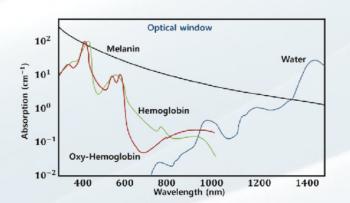

FEATURES

OPT / AFT Technology

Keep every pulse same and even, replace traditional sharp triangle output mode. It makes the penetration depth and energy more uniform, the procedure more comfortable, especially much safer for dark skin type.

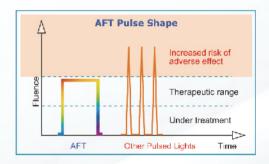

Two-step Filtering Technology


Make the wavelength much purer, more effective and painfree.

2 Working System. In-Motion & Stamp

In-Motion is up to 10 shots per second. Stamp for single shot.

2 Handpieces for Multi-function

German Xenon Lamp with the most scientific handle structure ensures the longer lifespan. Each handle is effective up to 1 million shots.

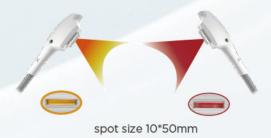

INTERFACE

Multifunctional and user friendly interface.







BEFORE & AFTER













SPECIFICATIONS

Screen Size	15"
Screen Size	· · ·
Handpieces number	2PCS
Wavelength	530nm - 750nm, 640nm - 750nm
Cooling System	Air+Water+Semiconductor
Dimension of machine	W45*D39*H30cm
Weight	56kg



ANDE SMART

Ande Smart represents a cutting-edge breakthrough in therapeutic innovation, seamlessly integrating advanced technology with intelligent design to deliver highly effective and personalized massage treatments. The artificial intelligence robot takes the place of human hands, offering services that are more precise, stable, safe, and efficient. This not only allows businesses to save on labor costs but also enables patients to save valuable time.







Body Shaping

Skin Tightening

Muscles Relaxing





WORKING THEORY

Equipped with an artificial intelligence probe, utilizing real-time monitoring and positioning technology, fusing cutting-edge radiofrequency energy, Ande Smart achieves automatic and precise control over high-standard temperature, output, pressure, and the automated generation of massage trajectories. Furthermore, its unique depth camera technology further enhances the accurate identification of treatment areas, enabling deeper and more precise penetration.



- Intelligent mechanical arm design breaks through traditional models, targeting every inch of skin and breaking the limits of traditional treatments.
- Deep tissue penetration and precise targeted therapy provide effective relief.
- Intelligent sensing for personalized care adapting to individual body contours. Live adjustment of massage intensity ensures a comforting experience.
- Allow for personalized treatments from soothing massages to deep repairs, help regulate sub-health conditions, and boost vitality.
- Incorporates ergonomic, skin-friendly design, and real-time collision detection to ensure a secure and efficient therapy experience.

BENEFITS

- Real-time detection, rapid response, and ensuring personal safety.
- Adjust massage intensity in real time to ensure individuals receive precisely customized pressure levels.
- The effect exceeds that of manual labor, providing deeper penetration, greater strength, and enhanced operational stability.
- Provide personalized, operator-independent services, enhancing the overall user experience.
- The heat therapy function regulates sub-health status, accelerates microcirculation, relieves fatigue, and effectively replenishes vitality.

SPECIFICATIONS

Size	18.5" Touched Screen	
Resolution	1920*1080	
Touch	Capacitive Touch	
Input Voltage Range	200-240VAC	
Frequency Range	50/60 HZ	
RF Energy	1-100W	
RF Frequency	1Mhz	
Machine Size	611*535*1254	
Machine Package Size	830*780*1640	
Machine + Package	150KG	
Machine Net Weight	87KG	

360°CRYOLIPOLYSIS COOL 360

Cryolipolysis, commonly referred to as fat freezing, is a nonsurgical fat reduction procedure that uses cold temperature to reduce localized fat deposits or bulges that do not respond to diet and exercise.



Accurate
Temperature Control



360° Surrounding Cooling



Non-invassive



WORKING THEORY

Cryolipolysis is a well-known slimming technology for permanent fat removal. It is using controlled low temperature to freeze fat and cause cell death of subcutaneous fat tissue, then the body will naturally eliminate the dead cells later via lymphatic system. Fat cells are more sensitive to cold temperatures than other types of surrounding body cells, so it has no effect on normal skin tissue.



Many of us have stubborn fat despite diet & exercise



Cryolipolysis technology uses controlled cooling to target & kill only these fat cells.



In the weeks to follow, your body naturally processes the fat and eliminates these dead cells.



Cryolipolysis procedure results are long-term, as treated fat cells are gone for good with no rebound.

APPLICABLE AREAS

Double chin, Thighs, Abdomen, Flanks, Bra fat, Back fat, Underneath the buttock, Upper arms.

INTERFACE

Intuitive and user-friendly interface for male and female.

Customized parameter for different body parts. Multi languages menu.






APPLICATOR

COOL 360 comes with 4 different size applicators for different body areas, and 9 cups for different fat thickness.

4 applicators can work seperately or simultaneously for different body areas or different patients.



The latest 360° surrounding cooling technology brings the uniform and even cooling of the whole area, which will be more effective than traditional two sides cooling technology.

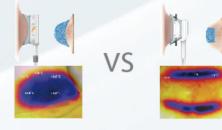

MUSCLE BUILDING MACHINE

EM Contouring can treat various parts of your body including abdomen, buttocks, arms, calves, and thighs to strengthen and improve muscle tone. The end results are an approximate 16% increase in muscle mass with a simultaneous average of 19% fat decrease.



NON-INVASIVE



OPERATOR INDEPENDENT



CONSUMABLES



DOWNTIME





.

Function

▼ Body contouring

✓ Muscle building



▼ Buttocks lifting

EM CONTOURING (Optional: RF technology)

EM MAX (Optional: RF technology)


WORKING THEORY

EM Contouring is based on HIFEM (High-Intensity Focused ElectroMagnetic) technology. A non-invasive procedure can induce 20,000 supramaximal muscle contractions in 30 minutes. These induced muscle contractions are very intense, therefore the muscle tissue is forced to adapt to such extreme conditions. The muscle tissue responds by deep remodeling of its inner structure, which results in muscle building and body contouring.



2 Different Size Applicators For **Different Body Parts**



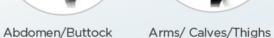













SCIENTIFIC MAGNETIC PULSES IN 3 WORKING STAGES

Warm-up and muscle training stage. Muscle training and fat cell apoptosis stage. Detoxification and relaxation stage.

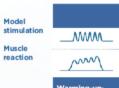



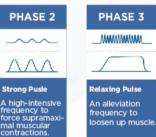

4 WORKING MODES

- Simple 1 For people who seldom exercise with less muscle and more fatness.
- Simple 2 For people who do regular exercise with more muscle and less fatness.
- Specialty 1 For professional practitioners to customize protocol A.
- Specialty 2 For professional practitioners to customize protocol B.

1060NM DIODE LASER SHAPE LASERSHAPE - A

Even with diet and exercise, stubborn fat seems impossible to lose. LaserShape's light-based body contouring technology targets and destroys fat in problem areas such as the belly and love handles, helping you achieve a slimmer appearance with no rebound.



25 Minutes Treatment



6-8 Weeks Up To 24% Permanent Fat Reduction



No Surgery. Non-Invasive




No Pain. No Downtime



WORKING THEORY

1060NM Diode Laser Shape has a high affinity for subcutaneous adipose tissue. The energy transmitted by the laser causes the molecules in the irradiated tissue to move, create heat, then raise the temperature between 42 °C and 47 °C. Around this temperature, the cell membrane structural integrity of the target adipocyte will be broken down, then leads the cells inactivation. During the next several months, the dead adipocytes can be eliminated by lymphatic metabolism, thereby achieving the effect of slimming and shaping.



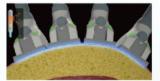

INTERFACE





The colors of A B C D will change smartly to show different working stage duing the whole procedure. Automatic Alarm System for safety.


APPLICATOR





4 Applicators (able to work at the same time) Spot Size: 40*80mm

Freely combine applicators to target areas you where want to get rid of stubborn fat.

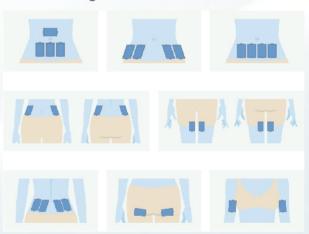

SPECIFICATIONS

15"
1060nm
W59*D61*H145cm
97kg

MULTISHAPE BODY SLIMMING SYSTEM VS++

Multishape is the evolution of the slimming and skin tightening technologies through the combination of the Powerful Cavitation, Multi-Polar RF, Vacuum and Cavitation. Using the Super-Pulse Technology and slidding motion from the Vacuum to achieve skin renewal, cellulite reduction and body reshaping management.

Function

M Body remodeling & Reshaping

✓ Cellulite & Fat reduction


▼ Facial lifting & Skin tightening

✓ Improve blood circulation

✓ Speed up body metabolism



UNIQUE WORKING PROCESS

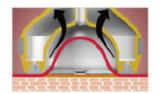



Suction

Release

15 pulses output per second to achieve maximum effect of physical activity.

Able to control the level of Suction and Release to customize different needs.

INTERFACE

Intuitive and user-friendly.





APPLICATORS

Using a diverse array of applicators, it is an easy and convenient way for whole body areas, face, arms, legs, abdomen, back and buttocks.









EFFECTS OF MULTI-POLAR RF

Increases intracellular diffusion of oxygen by heat.

Increases the microcirculation.

Improves the elasticity and firming of the skin.

EFFECTS OF VACUUM

Relaxes the tissue and massages the fascia Relieves swelling and greatly improves body and facial contour line.

EFFECTS OF ULTRASONIC

Fat promotes degradation of fat by reaching deep into the area of subcutaneous fat.



LARGE HANDLE RF & Vaccum & LED Light (Red)





SMALL HANDLE RF & Vaccum & LED Light (Blue)



ULTRASONIC CAVITATION HANDLE







MULTISHAPE BODY SLIMMING SYSTEM VS+

FOR BODY

Body Slimming & Firmness Cellulite & Fat reduction Waistline loss Body massage

FOR FACIAL

Firmer skin Wrinkles smoothed Firmed facial contours


FOR HEALTHCARE

Relieve fatigue Accelerate Metabolism Reduce pain on neck and low back



WORKING THEORY

VS+ combines 4 technologies of Vacuum, RF, Infrared Laser, Mechanical Roller to heat the fat cells, their surrounding tissues and the underlying dermal collagen fibers.

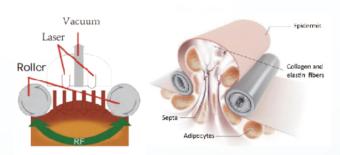

This non-invasive procedure can contour and slim the body by reducing cellulite, disruptng localized fat deposits, and tightening and firming skin.

INTERFACE

Intuitive and user-friendly. Able to save 5 groups parameters.





APPLICATORS FOR DIFFERENT **BODY PARTS**

HANDPIECE 1

Vacuum + Roller + RF + Infrared light For treatment on abdomen and back







HANDPIECE 2

Vacuum + Roller + RF + Infrared light For treatment on shoulders, arm ,abdomen, back, buttocks, thighs







HANDPIECE 3

Vacuum + RF + Infrared light For treatment on face and neck







SPECIFICATIONS

Screen Size	15"
Handpiece number	3pcs
Dimension of machine	W76*D51*H144cm
Weight	56kg

LEGENDA 448KHZ DEEP BEAUTY FOR BODY CARE

PHYSIO AESTHETICS

Visceral fat removal Fat and cellulite reduction Face & body contour Skin tightening & rejuvenation



REHABILITATION

Osteoarthritis Low back pain & Neck pain



SPORTS PHYSIOTHERAPY

Reduce pain Haematoma Muscle tears Sprains/Strains



Fractures



WORKING THEORY

This Legenda converts an electric current into a unique and stable frequency of 448KHZ to stimulate tissue repair. Extensive scientific research proves that 448 KHz is the optimal frequency to produce the best therapeutic outcomes today for Wellness and Aesthetics. It delivers the energy in the form of CAP and RES to target the different depth of tissue, achieving the effects of accelerating the recovery from injuries, reducing pain of musculoskeletal problems, improving pelvic health after child birth or or aged skin, and enabling aesthetic purposes of skin firming and body contouring.





INTERFACE

Intuitive and user-friendly. Automatic energy detection.





BEFORE & AFTER





2 WORKING MODES CAP & RES

CAPACITIVE MODE (CAP)

Suitable for superficial and vascularised tissue

The reactions caused by capacitive mode are mainly focused on the soft tissue adjacent to the electrode. Its design enables it to act as an external capacitor and to obtain a localised increase of electro-regenerative energy. The capacitive mode reduces the internal resistance of the tissues to the passing of current (as a result of the improved vascularisation).



RESISTIVE MODE (RES)

Suitable for thick, dense, fatty and fibrotic tissues

The reactions caused by the resistive mode act at a deeper level and focus on more resistive tissues, such as joints. ligaments and tendons. This mode in combination with manual techniques and kinetic therapy enables concentration of the effects of electro-regenerative energy on these structures.







SPECIFICATIONS

	Screen Size	8"
	Frequency	448KHz
	Dimension of machine	W49*D50*H37cm
	Weight	22kg

Take a sit and experience the breakthrough technology for incontinence.

30 Mins **= 11,000** Kegels



NO SURGERY



NO PILLS



NO PAIN



Function

✓ Strengthen the entire pelvic floor muscles


✓ Improved sexual performance

Reduced incontinence

REDUCED INCONTINENCE

EM Bella utilizes High-Intensity Focused Electromagnetic technology (HIFEM) to deliver thousands of supramaximal pelvic floor muscle contractions in a single session. These contractions re-educate the muscles of incontinent patients similar to the contractions you perform when doing a Kegel exercise. EM Bella is one of the easiest, most relaxing, and non-invasive procedures available to help reduce urinary leakage and dribbles, and free women and men from incontinence.

Myogram

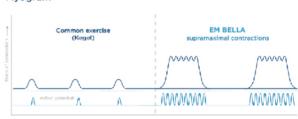

■ FFMALF

■ MALE

60% female are affected by UI between the ages of 45 and 55

15% males are affected by UI over age of 60



IMPROVED SEXUAL PERFORMANCE

Erectile dysfunction, or ED, a signififcant and common medical problem, is the inability to achieve or sustain an erection suitable for sexual intercourse.

35% of men aged 40-70 are suffered from ED 5% -10% of men below the age of 40 are suffered from ED

The EM BELLA is one of the most effective ways to combat ED by using electromagnetic pulses to stimulate the musculature of pelvic floor. These pulses help restore neuromuscular control, which can help prevent incontinence but also work as an exercise to reverse ED.



BEFORE

Pelvic floor muscles insufficiently support pelvic organs and affect bladder control.



PROCEDURE

EM BELLA effectively stimulates pelvic floor muscles with thousands of contractions per session



AFTER

Pelvic floor muscles sufficiently support pelvic organs and better bladder control

SPECIFICATIONS

Model	EM Bella	EM Bella Pro	
Screen Size	/	15"	
Energy	3 Levels	1% - 100%	
	adjustable	adjustable	
Voltage	110V - 220V	110V - 220V	

Q-SWITCH ND:YAG LASER

Function

✓ All colors tattoo removal

☑ Lip line & eyeliner removal

✓ Pigmented Lesions removal

Skin whitening & skin rejuvenation

✓ Carbon peel skin rejuvenation



WORKING THEORY

The Q-switched ND YAG laser takes specific wavelength light in high energy, which is absorbed by the pigment and shatters the pigment into particles, breaking them into very small fragments, some parts will consequently bounce out of the skin and the other parts will spilt even further into minute particles which eventually will be engulfed by the phagocytes and ultimately get eliminated by the lymphatic system.







Due to photoacousic effectpigment particles are selectively by phagocytes
 Fragmented pigment particles are carried away by phagocytes

INTERFACE

- 10"Display
- 5 Multiple optional functions
- Single pulse -1300mj
- Double pulses 1600mj





DUST-FREE ROOM

China largest Dust-free Room for production of optical laser.



LASER GENERATOR

Patented technology

The most scientific structure to achieve high power and stable energy output Easy Calibration



KOREA ORIGINAL LASER ARM

Korea imported 7-joint laser arm Light transmission loss less than 10%



2 Standard Handpiece







Skin Rejuvenation Handpiece

SPECIFICATIONS

10.4"
1064nm&532nm
1600mj
220V/110V

(SUPER) PICOSECOND LASER

Picosecond presents a new aera in deep pigments shattering, offering optimal results for the removal of colorful tattoos, pigmented lesions, wrinkles and acne scars with just few sessions and virtually no downtime.

Function

▼ Body tattoo

▼ Eyebrow tattoo

✓ Eyeliner tattoo

▼ Melasma

✓ Naevus

▼ Freckles

Age spot

✓ Coffee spot

Nevus of Ota


Skin toning


✓ PIH (Postinflammatory) hyperpigmentation)



WORKING THEORY

Picosecond mode of operation is based on delivering ultra-short picosecond pulses of energy into the skin tissue by the principle of light mechanical shock wave. Then the pigment is "shattered" via focused energy into ultra-fine granules, which are more likely to be absorbed by the body metabolism. Picosecond laser can minimize the side effects of thermal effect, and achieve the goal of solving almost all kinds of pigment spots, it is better than traditional nanosecond laser.

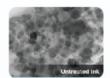

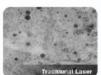



BEFORE & AFTER



INTERFACE

- Multiple optional functions
- Automatic protection system





KOREA ORIGINAL LASER ARM

Korea imported 7-joint laser arm Light transmission loss less than



2 STANDARD HANDPIECE





1064nm&532nm Handpiece

Honevcomb Handpiece

SPECIFICATIONS


	PS2014	PS10	
Screen size	10.4"	15"	
Wavelength	1064nm; 532nm	1064nm; 532nm	
Pulse width	750ps	450ps	
Frequency	1-10Hz adjustable	1-10Hz adjustable	
Spot size	2-10mm adjustable	2-10mm adjustable	

ULTRAPULSE FRACTIONAL CO2 LASER



WORKING THEORY

The Ultrapulse Fractional CO2 laser emits light with a wavelength of 10600nm, allowing the evaporation of microscopic parts of the skin, leading to heating of deep structures. Thermal damages strongly stimulate the synthesis of new collagen, its remodelling, the formation of new regenerated cells.

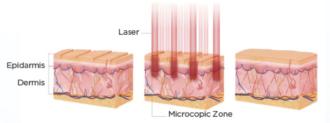

3 WORKING MODES WITH MULTI-FUNCTION



- Fractional Mode
- Collagen remolding
- Stretch marks removal
- Scars removal
- Shrink pores
- Skin tightening & resurfacing & rejuvenation




- CO2 Surgery Mode
- Warts removal
- Moles removal
- Skin tags removal



- 360°Mode
- Vaginal Tightening
- Vaginal Dryness
- Vulva Rejuvenation

7 SCANNING PATTERNS

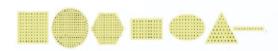

3 SCANNING MODELS







Sequence

Maxium Distance

Random Scan

INTERFACE







SPECIFICATIONS

CO2 Laser Source	RF Mental Laser Tube
Screen size	10"
Wavelength	10600nm
Dimension of machine	W44*D44*H112.5cm
Weight	56kg

GOLDEN AGE MICRONEEDLING RF







LIFTING



ELASTICITY



Collagen remolding

Stretch marks removal

Scars removal

Shrink pores

Skin tightening & resurfacing &

rejuvenation



WORKING THEORY

Microneedling RF (MNRF) creates tiny skin perforations and delivers thermal energy to underlying tissues to cause dermal coagulation, collagen remodeling and neoelastogenesis via a wound-healing response.



INTERFACE

Intuitive and user-friendly interface. Customized for different body parts. Multi languages menu.





HANDPIECE

Ergonomic handpiece design Insulated needles







PRECISE DEPTH CONTROL

0.5-3.5mm ajustable depth with 0.5mm step.

BEFORE & AFTER





SPECIFICATIONS

Screen size	10"
Working mode	Auto/ Manual
Dimension of machine	W44*D44*H125cm
Weight	42kg

ME-TIME

ME Time - Simultaneously tone facial skin & muscles.



NON-INVASIVE



Quick Results



Hands-Free



NO DOWNTIME

Function

▼ Facial Muscle Enhancement


✓ Wrinkles Reduction


✓ Skin Lifting



WORKING THEORY

ME Time works on the framework of the face from the inside out and targets facial layers with cutting-edge tech. It restores and elevates the support of facial tissues by selectively stimulating muscles and increasing the density and quality of muscle structure, enhancing the facial structure, and simultaneously remodels and smooths the skin by increasing collagen and elastin fibers, which results in a reduction of wrinkles and fine lines while tightening the skin.

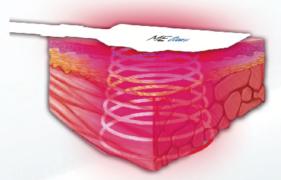

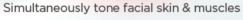



FEATURE

- One course for overall improvement of the structure of the face and skin quality.
- Gentle and pain-free non-invasive facial treatment requires no scalpels, no needles, no fillers.
- Safe for most people, including those with sensitive skin.
- Unique technologies to target all facial layers.
- No downtime.

BEFORE & AFTER









360 ULTRA HIFU HIGH INTENSITY FOCUSED ULTRASOUND

360 Ultra HIFU, a non-invasive treatment, goes deeper to stimulate production of your body's own collagen and elastin, for natural-looking results without surgery or downtime.











FOR FACIAL

Facelifting Skin tightening Wrinkles removal Collagen regeneration Improve skin elasticity



FOR BODY

Body slimming Fat loss Buttocks lifting

HANDLE













INTERFACE







Automatic detection when installing the cartridge of different depth.

Customized parameters for different body areas. Able to check remaining shots and shots of a single session.

FEATURES

Accurate depth control. Long lasting results after a single session. Non-invasive. No downtime.



BEFORE & AFTER





SPECIFICATIONS

Screen Size	15"
Handle number	3pcs
Dimension of machine	W40*D45*H60cm
Weight	30kg

Function

▼ Facelifting

✓ V-Shaped Face

✓ Skin Tightening

✓ Wrinkles Removal

▼ Eye Bags Removal

✓ Dark circles Removal

▼ Body Slimming



WORKING THEORY

20MHz Cold Capacitive-coupled Radio Frequency (CRF) Output "+" "-"electricions to generate movement, whirling, friction, collision with the "+" "-" electricions in human body, causing 65-70°C thermal energy to heat skin, causes collagen to contract and re-form, and triggers the production of new collagen and elastin firming, smoothing, lifting the skin and reducing the appearance of lines.







UNIQUE COOLING & HEATING WORKING METHOD

Contact cooling on the skin surface and heating under the deep layer, which provides a comfortable feeling much like a warm massage. Cold CRF will provide you with a long-lasting and immediate results without any pain and downtime.



INTERFACE

Intuitive and user-friendly interface. Customized for different body parts. Multi languages menu.





HANDPIECE

Ergonomic handpiece design to prevent hand fatigue for long time operation.

3 size tips (18mm,30mm,50mm) for different body areas.





SPECIFICATIONS

10"
20MHz
1PCS
W49*D49.5*H108cm
39kg

ULTRA TONE MULTIFUNCTIONAL 2 IN 1 SYSTEM

FACE & BODY

- Facelifting
- V-Shaped Face
- Skin Tightening
- Wrinkles Removal
- Eye Bags Removal
- Dark circles Removal
- Body Slimming

VRL

- Vaginal Tightening
- Vaginal Dryness



WORKING THEORY

20MHz Cold Capactive-coupled Radio Frequency (CRF) Output "+" "-"electricions to generate movement, whirling, friction, collision with the "+" "-" electricions in human body, causing 65-70 °C thermal energy to heat skin, causes collagen to contract and re-form, and triggers the production of new collagen and elastin firming, smoothing, lifting the skin and reducing the appearance of lines.







UNIQUE COOLING & HEATING WORKING METHOD

Contact cooling on the skin surface and heating under the deep layer, which provides a comfortable feeling much like a warm massage. Cold CRF treatment will provide you a long-lasting and immediate results without any pain and downtime.







HANDPIECE

- ① Handpiece with 3 sizes tips for cantus, facial and body.
- ② For vaginal function.



BEOFRE & AFTER





SPECIFICATIONS

Screen Size	10"
Frequency	20MHz
Handpiece number	2PCS
Dimension of machine	W45*D53*H36cm
Weight	28kg

Function

✓ Multi-functional

✓ Cost-effective

✓ Essential machine for Salons/Spas



WORKING THEORY



Cooling & Heating Probe

Improve skin redness and sensitivity



Ultrasound

Ultrasound can activate fibroblasts to improve elasticity and appearance.



ION Lifting

It uses micro current to stimulate underlying skin tissues and muscles for a more youthful appearance.



Multipolar RF

The optimum RF energy stimulates collagen and elastin to tighten skin, also improve blood circulation and lymphatic system.



Aqual Peel

It adopts vacuum suction and eddy technology, the spiral design tip can easily dislodge and remove impurities, and simultaneously hydrate skin by using solution as the medium.



Oxygen Spayer

This spraying process can deliver 98% pure oxygen infused with serum (water-soluble products) in super- nanoparticles to the skin tissue, getting better absorption and hydration.



BORE Handle

This treatment can create BOHR effect, combining exfoliation, nutrient infusion, and oxygenation from within.



STEP 1 EXFOLIATION

It removes the upper layer of dead skin cells, similar to a gentle microdermabrasion.



STEP 2 INFUSION

essential nutrients and antioxidants are infused deeply into the skin, skin permeability is increased.



STEP 3 OXYGENATION

It produces carbon dioxide bubbles that gently burst on the skin's surface, which then stimulates oxygen-rich blood flow to the area and greater absorption of the actives.

BEFORE & AFTER





SPECIFICATIONS

Screen Size	10"
Handpieces number	7PCS
Dimension of machine	W45*D39*H30cm
Weight	31kg


QUEEN LIGHT BIO LIGHT SKIN REPAIRING SYSTEM

QUEEN LIGHT (also called LLLT or LPLT) has been a FDA-cleared skin care project since 2003. It can help us resolve complicated skin problems such as acne, acne scars, pigmentation, dull skin, sensitive skin and more.



WORKING THEORY

QUEEN LIGHT uses excellent and strong high-brightness special near-infrared SMD LED light with various wavelengths. The light emitted is proven to be beneficial in cellular activity and repair, dramatically revitalizing skin without any downtime or discomfort.

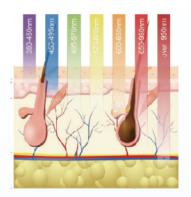

7 COLORS FOR MULTIFUNCTIONS



Color	Wavelength	Function
Red Light	622nm	Anti-Aging
Blue Light	462nm	Acne
Green Light	522nm	Against Hyperpigmentation
Yellow Light	432nm	Against Sunburn and Other Redness
Purple Light	400nm	Eliminate Acne and Redness
Cyan Light	490nm	Double Anti-acne and Anti-hyperpigmentation
White Light	400-700nm	Accelerate Tissue Metabolism

INTERFACE

User-friendly operation interface Multiple languages menu









SPECIFICATIONS

Light source	Bio Light
Working mode	Continous/ Pulse
Colors	7 colors
Dimension of machine	W44*D44*H125cm
Weight	42kg

HYDRA MAGIC MESOTHERAPY INJECTION SYSTEM












WHAT IS HYDRA MAGIC?

Hydra Magic is an injection-supportive device to be designed to help the users to operate mesotherapy easily and accurately with stimulation of mesoderm by using syringe with drugs.

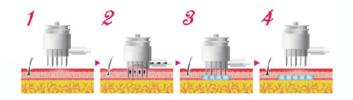

DRUGS FORMULATION

HA + POLYMER HYALURONIC ACID (HA)

Reduce fine lines and wrinkles, moisture skin, hydrating, increase the elasticity of the skin.

HA + VITAMIN C

The best anti-aging ingredients, and the key to maintaining a smooth, even, and glowy complexion.

HA + PLACENTA SERUM

Be used to strengthen and re-build the skin. Ideal for post-acne scarring and after peels or any form of epidermal exfoliations.

HA + PDRN

PDRN is a naturally-occurring building block in your DNA that helps to repair damaged skin cells and stimulate collagen.

INTERFACE





FEATURES

AUTOMATIC INJECTION SYSTEM

Auto-suction, Auto-injection, Auto-release

ACCURATE CONTROL OF DOSAGE

The minimal dosage can be accurate to 0.0083cc

LIGHT ERGONOMIC HANDLE DESIGN

Comfortable hold, easy operation

SMALL AND EASY TO CARRY

Save space and convenient for door service

KOREA MADE NEEDLES

High hardness, Less painful, No drug leaking





BEFORE & AFTER




SPECIFICATIONS

Pressure	~620mmHg
Dimension of machine	W26*D32*H15cm
Weight	2.5kg

VIEWSKIN 4D SKIN ANALYZER

3D Skin Analyzer can provide a complete analysis of all skin problems, not only the problems of epidermis, but also dermis problems. Clients can see their skin problems in a scientific method and the staff will sell rejuvenation and skin care solutions more effectively with clear recommendations combined with the visual impact of our 3D Skin Analyzer.

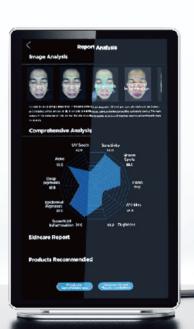



WORKING THEORY

3 Light Sources (Standard light, Cross Polarized light, Ultraviolet light) with a 35 million pixels industrial camera.

















To ensure the accuracy of analysis.



- Pores
- Wrinkle
- Acne
- Sensitive
- Roughness
- Deep pigmentation
- Superficial inflammation
- Epidermal Pigmentation














FEATURE

operation, easy operation, intuitive, sewing patterns at a glance.

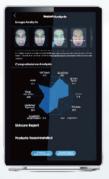

SMART CLIENTS DATA SYSTEM

All clients' data can be archived in the system, which helps to record and manages all your clients' datas effectively.

SPECIFICATIONS

Product name	ViewSkin
Pixels	35 million
Light Sources	Standard light
	Cross Polarized light
	Ultraviolet light
System Language	English, Chinese



VIEWSKIN MEDICAL SKIN ANALYZER



COMPANY PHILOSOPHY

MISSION

Lead the development of the global beauty industry with honesty and morality.



VISION

To be a leading enterprise in the global beauty and health industry.



VALUE

Patriotism Integrity Morality Innovation Responsibility Progress



CONCEPT

Do only valuable things for customers.