

HemoCue® HbA1c 501 System



Designed to Work for How You Work

HemoCue HbA1c 501 is a fully automated point-of-care system that's not only easy to use, but every element – from its size to the storing and handling of cartridges – is thought through to be ideally suited for how you work.

Rigorously tested according to HemoCue principles, and with automatic checks and calibrations built into the system, you can count on getting the right answer at the right time – every time.

So you won't miss a beat in caring for your patients.

Accuracy Starts With Us



To learn more about HemoCue® HbA1c 501, please scan the QR-code with your smartphone or visit hemocue.com

Makes Your Whole Process Easier and More Efficient

- ▶ Only a few simple steps with clear instructions on screen
- ▶ All-in-one, individually packaged cartridge kits
- ▶ Convenient room-temperature storage up to 18 months

Ensures You'll Never Have to Doubt an Answer at the Point of Care



- ▶ Interference-free method unaffected by Hb variants
- ▶ IFCC and NGSP certified
- ▶ Precise, fixed calibration by factory

Provides a Safe Workflow for Excellent Patient Care


- ▶ Automatic reading of calibration data for each test cartridge lot
- ▶ Barcode reader for patient and operator ID available
- ▶ Connectable

HemoCue® HbA1c 501 System


Principle	Fully automated boronate affinity assay for determination of HbA1c percentage in whole blood
Calibration	Factory calibrated and traceable to IFCC and NGSP/DCCT
Sample Material	Capillary or venous whole blood
Measurement Range	20–130 mmol/mol (IFCC), 4.0–14.0 % (NGSP)
Results	In 5 minutes
Sample Volume	4 µL
Dimensions	198 × 217 × 136 mm (7.80 × 8.54 × 5.35 inches)
Weight	1,600 g (3.53 pounds)
Storage Temp.	Analyzer: 10–35 °C (50–95 °F) Test Cartridge: unopened 2–32 °C (36–90 °F)
Operating Temp.	17–32 °C (63–90 °F)
Power	9 V DC / 1.5 A
Interface	Printer, PC and Barcode Scanner
Quality Control	Built in "selftest", Daily Check Cartridge – to avoid errors during analysis, Monthly Check Cartridge - to ensure analyzer performance. The system can be verified using external quality control material

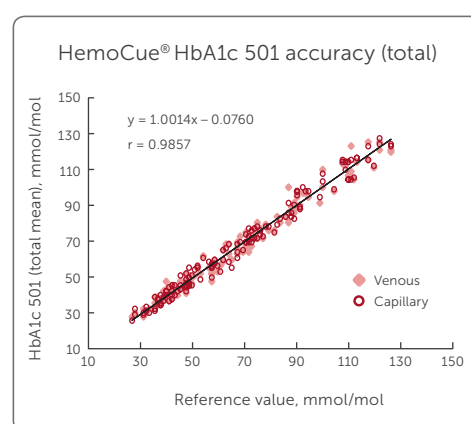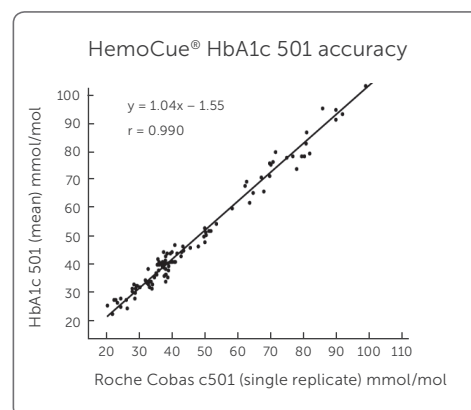

A few simple steps

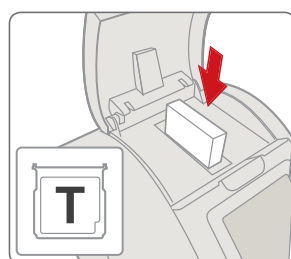


1 Insert the cartridge into the cartridge compartment.

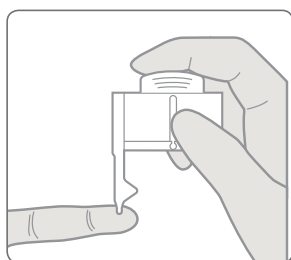

2 Prepare and apply the reagent pack to the specimen. The correct volume is drawn in by capillary action.



3 Insert the reagent pack into the cartridge. (Do not wipe off excess blood outside the sampling area.)

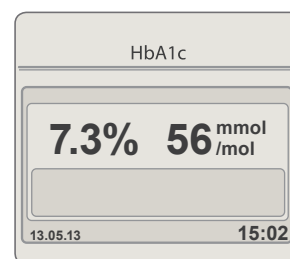

4 The lab-quality result is displayed automatically.

HemoCue AB | PO Box 1204 | SE-262 23 Ängelholm | Sweden
Phone: +46 77 570 02 10 | Fax: +46 77 570 02 12 | info@hemocue.se | hemocue.com
© 2015 HemoCue AB | HemoCue - part of the Radiometer Group

HemoCue has been a leader in point-of-care medical diagnostics for over 30 years. We specialize in giving healthcare providers lab-quality accuracy with results comparable to that of a clinical lab.

