



JR·BEAUTY


MAKE BEAUTIES MORE BEAUTIFUL



JR·BEAUTY

<http://www.jr-beauty.com>

PERMANENT Hair Removal

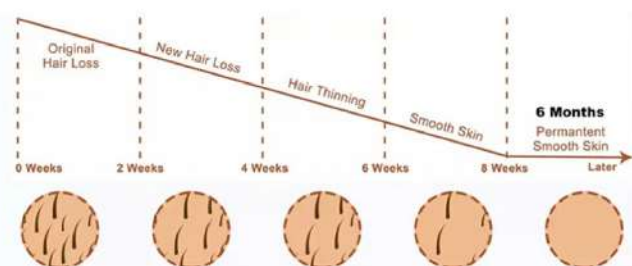

Product name	808+SHR+IPL+YAG 4IN1 Multifunctional Machine
808 Handle Laser bar	USA Coherent imported, 10 bar, 12 bars, 16 bars
808 Wavelength	755nm 808nm 1064nm
808 Laser handle shots	100 million shots times
808 Handle Laser Power	1000W/1200W/1600W optional
808 Spot size	13*16mm /13*22mm/13*30 optional
808 Frequency	1-10Hz
808 Pulse width	0-300ms
IPL Shots	UK imported lamp with 300,000 shots
YAG Shots	USA imported lamp with 5,000,000 shots

4IN1 MULTIFUNCTIONAL MACHINE


808+OPT+IPL+YAG

Gold Standard working energy

Super short pulse duration

Super long life span

Super high power

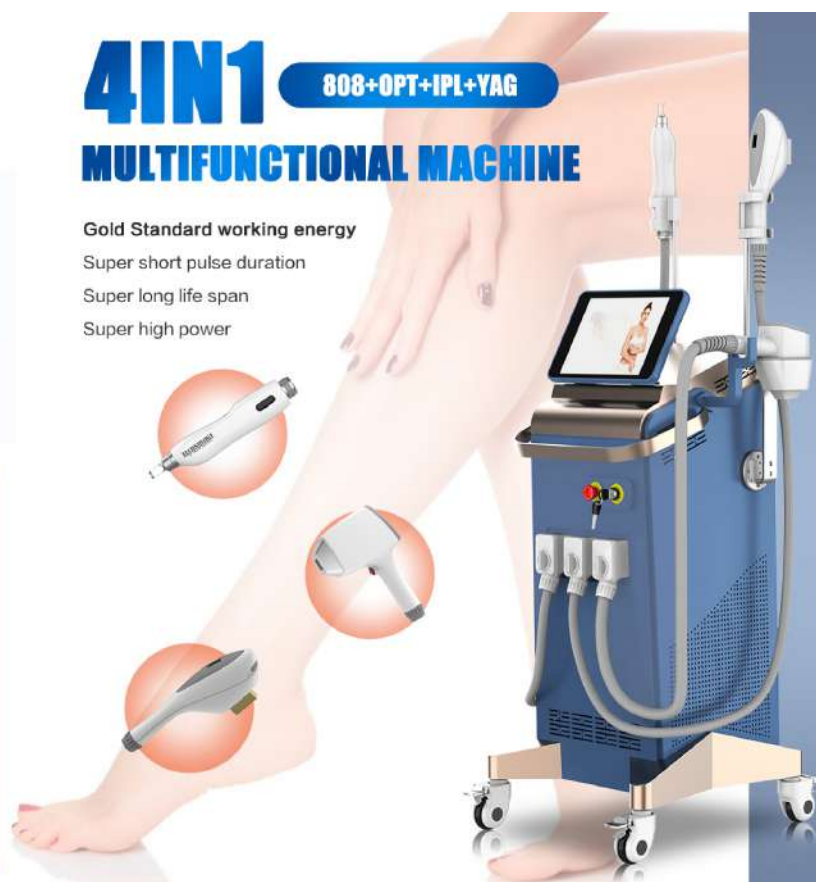



System with **16** different languages meet all customer demands

ADVANTAGES


- System chip can store 10,000+ customer information
- Thermoelectric cooler, better cooling system can make machine work continuously for 12 hours
- Double cooling plates can achieve real freezing point depilation
- System with 16 different languages meet all customer demands
- Water and electricity separation design, safe operation and zero failure rate
- Imported water pump nearly no noise



808nm DIODE LASER HAIR REMOVAL



Laser Type	Diode Laser
Laser Power Supply	3000w
Laser Stacks Power	1000w/1200w/1600w
Wavelength	808nm/1064nm/755nm
Spot Size	13*16mm13*22mm/13*30mm
Display	10.4" True Color LCD Touch Screen
Frequency	1-10Hz
Electrical Requirements	90-130V 50-60HZ or 200-260V.50HZ
Machine Size	520*380*370mm(L*W*H)
Net Weight	30kgs

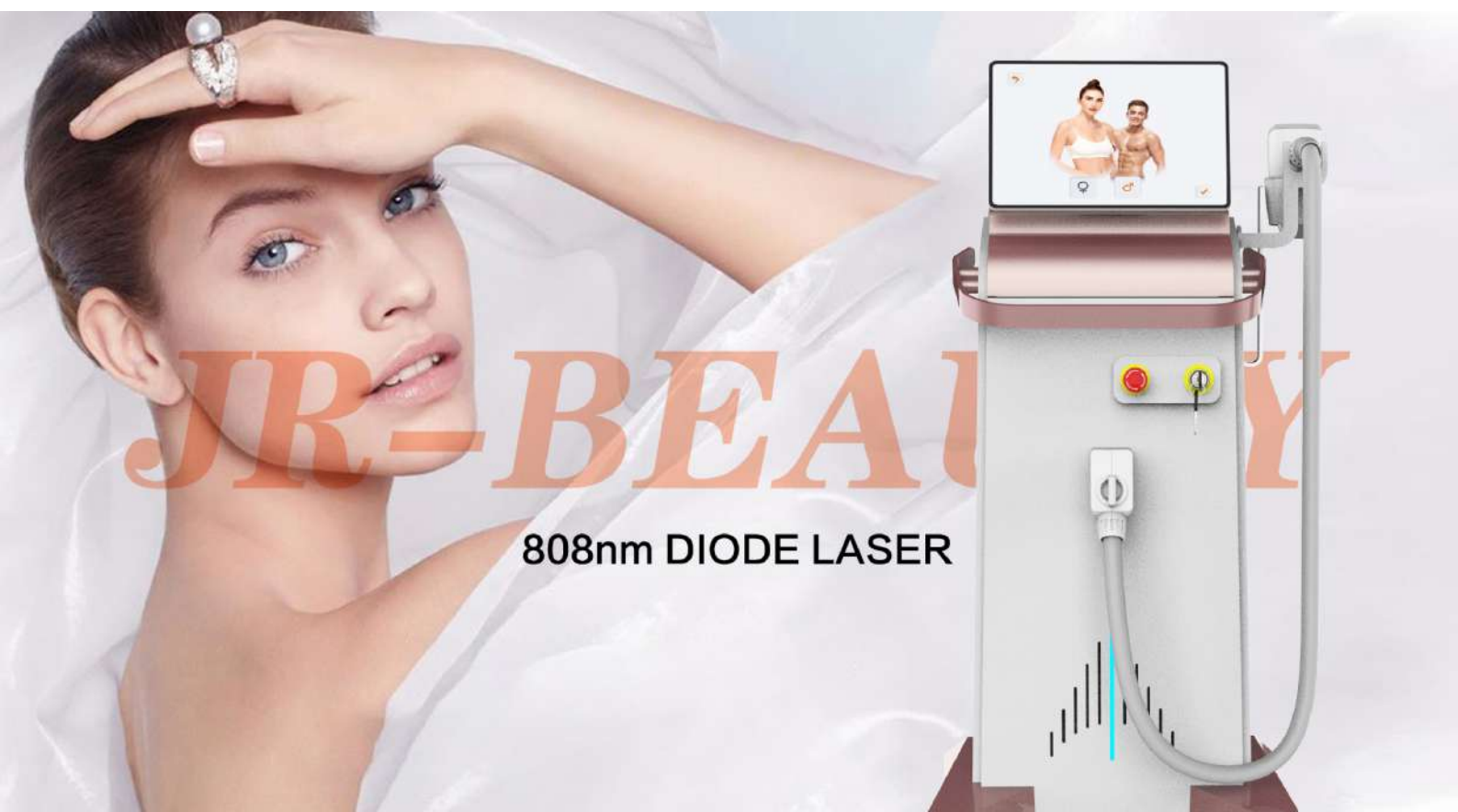

808nm DIODE LASER





Specification

Laser Type	Diode Laser
Laser Power Supply	4000w
Laser Stacks Power	1600w
Wavelength	808nm/1064nm/755nm
Spot Size	13*30mm
Display	10.4" True Color LCD Touch Screen
Manipulation language	System Language: Default English(Can be customized)
Electrical Requirements	90-130V.50/60HZ or 200-260V50HZ



Ice Cooling Effect



3 Highly Effective Wavelengths



Separation of water and electricity



All Skin Color



Permanent Hair Removal



Multi-size spot



Safe Painless



Mode selectable men and women



2-years Warranty



1 at your initial treatment the hair should be visible on the surface of the skin.



2 As the laser is directed onto the hair heat energy is drawn into the melanin, which in turn heats the rest of the hair shaft.



3 The heat either weakens or completely destroys the hair follicle. As a result the follicle detaches from the hair shaft.



4 The hair sheds and, without the follicle, hair cannot regrow in the area.

Laser Type	Diode Laser
Laser Power Supply	3000w
Laser Stacks Power	1000w/1200w/1600w
Wavelength	808nm/1064nm/755nm
Spot Size	13*16nm/13*22nm/13*30nm
Display	10.4" True Color LCD Touch Screen
Manipulation language	System Language: Default English(Can be customized)
Electrical Requirements	90-130V.50/60HZ or 200-260V50HZ

GERMANY LASER CHIP

- High energy density
- At least 100 millions shots
- Stable and longifespan
- 2 years warranty



1200W/1600W/2000W 808 DIODE LASER HAIR REMOVAL

Ice Cooling Effect

All Skin Color

Safe Painless

3 Highly Effective Wavelengths

Permanent Hair Removal

Mode selectable men and women

13*30mm

13*22mm

13*16mm

THE INTELLIGENT SOLUTION FOR HAIR REMOVAL



Double handles New **IPL** Laser Hair Removal Machine



Hair removal



Pigment removal



Acne removal



Vascular lesion




Skin rejuvenation



Skin tightening



IPL ELIGHT SKIN TREATMENT

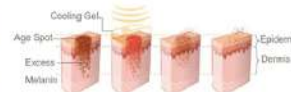

IPL technology is a type of Intense Pulsed Light that generates intense bursts of red or blue light. These heat the dermal tissue without harming or removing the outer epidermal layers of the skin. The body's natural processes then work to break down and remove these marks. The appearance of any skin imperfection is then visibly reduced.



Operation Interface	10.4 inch TFT Screen
Wavelength	610-950nm(HR)810nm/430nm/640nm/590 nm
Working Mode	Normal IPL filter working&3hit Working
Wavelength	15X50nm;15X50nm(option for both handles)
Spot Size	13x15mm;13x22mm/13x30mm
Energy	1-50J/cm2(SR)1-26J/cm2 (SHR & SHR)
Pulse Width	1-9.9ms
Electrical Requirements	90-130V 50-60HZ or 200-260V 50-60HZ
Pulse Amount	1-6
Handle	1X530-1200nm(SR)1X610-950nm(HR)

STANDARD SIZE



480nm



530nm



640nm

OPTIONAL SIZE



560nm



590nm



690nm



750nm

Super Hair Removal "In-motion" Hair Removal

Exposure to high pulse repetition rates at low energy levels gradually heat the hair follicle



Before

The melanin absorbs the light and gradually heats up

After



SUPER DPL IPL FILTERS LASER HAIR REMOVAL MACHINE



480nm



560nm



530nm



590nm




690nm



640nm



750nm

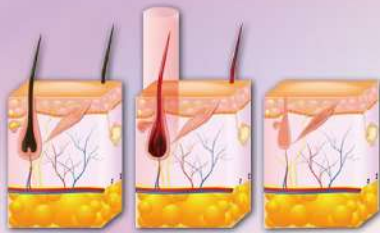

Before

The melanin ab-
sorbs the light and
gradually heats up

After

SPECIFICATION

Operation Interface	10.4 inch TFT Screen
Wavelength	610-950nm(HR)filter(430nm/640nm/590 nm)
Working Mode	Normal IPL filter working&SHR Working
Wavelength	15X50mm,15X50mm(option for both handles)
Spot Size	13*16mm 13*22mm/13*30mm
Energy	1-50J/cm2(SR)1-25J/cm2 (SHR & SSR)
Pulse Width	1-9.9ms
Electrical Requirements	90-130V 50/60HZ or 200-260V 50HZ
Pulse Amount	1-6
Handle	1X530-1200nm(SR)1X610-950nm(HR)





PICOTECH ND YAG LASER MACHINE

Advantage

- Less treatment times
- More effect on all pigmentation
- Avoid pigmentation risk
- Fractional head: No scab


LASER TYPE	ND: yag laser
Machine Power	1500W
Laser Power	800W
Output System	Handle Tool Output
Wave Length	1064nm&532nm&755nm
Pulse Width	532nm:4ns/755nm:4ns/1064nm:4.5ns
Energy	100-2000mj
Frequency	1-10hz
Cooling System	Wind Cooling+Water Cooling
Voltage	AC220V±10% 50Hz/AC110V±10% 60Hz

ND: YAG LASER Hand pieces

Hand pieces	Silver Hand piece Go With machine	Golden Hand piece GoWith machine
	1064nm	532nm
	Remove and lighten speckles, age- spots, sun- burn spots and tattoo in black, blue colors	Remove and lighten red andbrown colors tattoo

Tattoo Removal

Permanent Makeup and All Colors Tattoo Removal

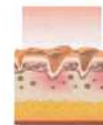

Pigmented Lesions Treatment

Nevus of Ota, birthmark, Sun Spot,Coffee Spot, Age Spot, Freckle



Skin rejuvenation

Oil control, Skin Whitening, Enlarged Pore Treat-ment, Carbon Laser Skin Rejuvenation



Smart water temperature:
Effectively protect the handle from damage



Unique shape design:
Products are more beautiful and fashionable



High spot uniformity:
Consistent treatment effect, safer




Short pulse width:
Quickly and effectively crush melanin



PRODUCT NAME	Q Switch ND YAG Laser
POWER	2000W
WAVELENGTH	1064nm/532nm(standard):755nm(optional)
FREQUENCY	1-10HZ
PULSE WIDTH	500-900ps
SPOT SIZE	2-10mm
DELIVERY SYSTEM	7 joints articulated arm
COOLING SYSTEM	Water + Wind + Radiator

Picosecond Technology

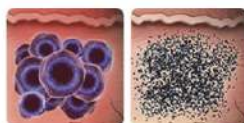


Picosecond takes advantage of pressure wave technology to shatter the target ink into tiny particles that are easily eliminated by the body.

Nanosecond Technology



Traditional nanosecond lasers predominately rely on photothermal action, delivering heat to the pigment and surrounding tissue.

Picosecond laser

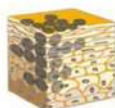

A Pigment hides in the dermis of the skin

Traditional laser




B Laser selectively acts on the pigment




1.Strongest energy:

Two benefits for stronger energy:a.expel the pigment faster and more thoroughly;b.with the decaying of the equipment energy, the stronger the energy is, the longer life of the equipment will be.

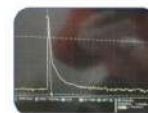

2.Shortest Pulse Width (500ps-900ps available):

Pulse width refers to the time that the laser applied to the skin, the shorter the time is, the less damage the laser leaves to the skin. It is also the reason of the popularity of the picolaser.

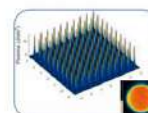

3.True " flat cap" output mode:


Uniform distribution of energy, each spot receives the same amount of energy. No energystrong point or weak point, which ensure energy being effectively absorbed by the target color, while prevent the normal skin from being damaged.



Top flat hat laser beam
Big and Perfect Round Spot



C Pigment blasts while absorbing the energy

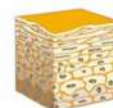


D Pigment is engulfed by macrophage and excreted

Picosecond laser tattoo removal machine


7 SCAN SHAPES

Adjustable Size And Density

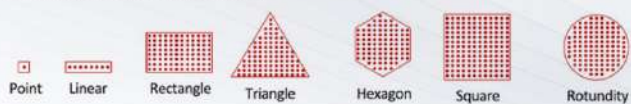

- Skin Resurfacing
- Skin Rejuvenation
- Stretch Marks
- Scar Remodeling
- Skin Sagging

Fractional Laser Skin Resurfacing

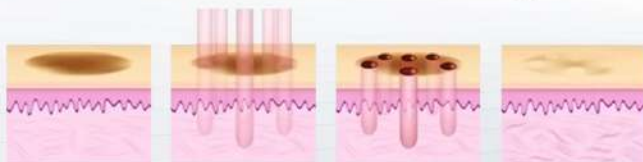

Vaginal Tightening



Laser irradiation
angle can be adjusted

The area of the laser
is compact edge contraction

Laser irradiation can be
rotated 360 degrees

Vaginal walls shrink
and tightening

Fractional *CO2* Skin Tightening Laser Beauty Machine



Removal of acne,
acne pits and acne marks



Treating old and
new scars



Vaginal treatment,
Tightening, Whitening



Treatment of pigmented
lesions such as freckles
and ostrom mirades




Anti aging, firm skin;
Remove fine lines from
the eyes and crow's feet



Eliminating pregnancy
induced wrinkles



Ultra Pulse Head mode

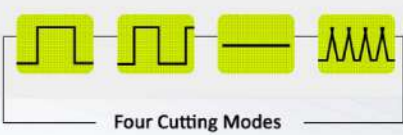

- Skin mole removal
- Flat wart removal
- Skin neoplasm removal

Vaginal Whitening

Vaginal Tightening

Normal mode

- Vaginal Tightening
- Vaginal Dryness
- Vaginal whitening
- Improve labial softness
- Post-delivery Rehabilitation

Output Power	1-60W (adjustable)
Input Power	450VA
Laser Type	Sealed-off CO2 laser
Cooling system	Inside-circle coolant
Divergence angle	4 mrad
Pulse width	0.1ms-10ms
Aiming beam	Red diode laser (650nm, 2mW)

- Remove freckles and Remove scars
- Removal of acne and acne marks
- Anti-aging: skin lifting and tightening
- Removing stretch of pregnancy
- Remove fine lines at the corners of the eyes and crow's feet
- Vaginal tightening, whitening





PRINCIPLE OF TREATMENT



FLAT HANDLE

- Abdomen
- Buttocks
- Thighs
- Back

CURVED HANDLE

- Arms
- Calves



Item Name	EMSlim Neo
Power	3000W
Voltage	100V-230V 50/60HZ
Pulse Width	300us
Intensity	13Tesla
Frequency	5hz-100hz
Muscle contraction	>20,000 times (with 30min)
Work Mode	UK imported lamp with 300,000 shots
IPL Opt size	Auto/ Manual



Treatment area: **ABDOMEN**



Treatment area: **ARM**



Treatment area: **BUTTOCKS**



Treatment area: **THIGH**

FOR THE BUTTOCK



FOR THE ABDOMEN

